Summary of Safety and Clinical Performance

Angio-Seal™ VIP Vascular Closure Device

Manufacturer	Terumo Medical Corporation 265 Davidson Ave, Suite 320 Somerset, NJ 08873 USA
Manufacturer's SRN	US-MF-000019594
Basic UDI-DI	38970AS6S
EMDN number	C90010302-Collagen-Based Haemostasis Systems
Device Class	Class III
Year Device was First CE Marked	1999
EU Representative	Terumo Europe, N.V. Interleuvenlaan 40, 3001 Leuven Belgium SRN: BE-AR-000001433
Notified Body	NSAI CE 0050

Table of Contents

1.0	Summary for Patients	4
1.1	Intended Use of the Device	5
1.1.1	Intended Purpose	5
1.1.2	Indications and Intended Patient Groups	5
1.1.3	Contraindications	6
1.2	Device Description	6
1.2.1	Description of the Device	6
1.2.2	Raw Materials & Components	7
1.2.3	Accessories Intended to Be Used in Combination with the Device	8
1.3	Risks and Warnings	8
1.3.1	Risks	8
1.3.2	Warnings and Precautions	9
1.3.3	Field Safety Corrective Actions	9
1.4	Clinical Evaluation and Post-Market Clinical Follow-Up	9
1.4.1	Summary of Data from Animal Testing	9
1.4.2	Summary of Data from Pre-Market and Post-Market Clinical Investigations	10
1.4.3	Clinical Background of the Device	44
1.4.4	Clinical Evidence Supporting CE Marking	44
1.4.5	Safety and Performance	44
1.5	Possible Therapeutic Alternatives	45
1.6	Suggested User Profile and Training	45
2.0	Bibliography	46
3.0	Revision history	48

This Summary of Safety and Clinical Performance (SSCP) is intended to provide public access to an updated summary of the main aspects of the safety and clinical performance of the Angio-Seal™ VIP Vascular Closure Device (VCD).

A summary of the safety and clinical performance of the device intended for patients is provided in **Section 1.0**.

1.0 SUMMARY FOR PATIENTS

The information presented below is intended for patients or lay persons. A more extensive summary of device safety and clinical performance prepared for healthcare professionals is found in **Section** Error! Reference source not found. of this document.

This SSCP is not intended to give general advice on the treatment of a medical condition. Please contact your healthcare professional if you have questions about your medical condition or about the use of the device in your situation.

This SSCP is not intended to replace an Implant Card or the Instructions for Use to provide information on the safe use of the device.

Manufacturer	Terumo Medical Corporation 265 Davidson Ave, Suite 320 Somerset, NJ 08873 USA
Manufacturer's SRN	US-MF-000019594
Basic UDI-DI	38970AS6S
EMDN number	C90010302-Collagen-Based Haemostasis Systems
Device Class	Class III
Year Device was First CE Marked	1999
EU Representative	Terumo Europe, N.V. Interleuvenlaan 40, 3001 Leuven Belgium SRN: BE-AR-000001433
Notified Body	NSAI CE 0050

1.1 Intended Use of the Device

1.1.1 Intended Purpose

The Angio-Seal™ VIP Vascular Closure Device is a device that is used to stop bleeding (hemostasis) from a hole in an artery in your groin (common femoral arterial puncture).

1.1.2 Indications and Intended Patient Groups

Indications for Use

The Angio-Seal™ VIP Vascular Closure Device is approved for use in closing a hole in the common femoral artery and help patients to walk around (ambulate) sooner after the procedure.

Target Patient Population

This device is meant for patients who have a hole in the common femoral artery and who are over 18 years old.

Angio-Seal™ VIP Vascular Closure Device has not been tested in the following patient groups:

- Patients who have known allergies to beef products.
- Patients with known disorder where your body's immune system mistakenly attacks your body (autoimmune disease).
- Patients getting blood thinning medication (therapeutic thrombolysis).
- Patients in which the hole was made through a device in the artery (stent or vascular graft).
- Patients with high blood pressure (uncontrolled hypertension, > 180 mm Hg systolic).
- Patients with a high risk of bleeding. Patients with a known bleeding disorder (thrombocytopenia, thrombasthenia, von Willebrand's disease) or the blood has low iron (anemia).
- Patients under 18 years old or with small arteries.
- Patients who are having a baby or breast feeding (lactating).
- Patients with blood pressure problems (Ankle-Brachial Index (ABI) of <0.9).
- Patients with blood sugar problems (uncontrolled diabetes).
- Patients with bleeding (hematoma, arteriovenous fistula, or pseudoaneurysm) at the procedure site.

- Patients experiencing shock (hemodynamic instability requiring intravenous medications or mechanical support).
- Patients with kidney or liver problems.
- Patients with other health problems (co-morbidity).
- Patients from a specific race or background (ethnicity).
- Patients who have had a heart attack (myocardial infarction) in the last 72 hours.

1.1.3 Contraindications

None.

1.2 Device Description

1.2.1 Description of the Device

The Angio-Seal™ VIP Vascular Closure Device is a medical device used to close a hole in the common femoral artery. It does this by making a sandwich like seal using plastic and animal tissue pieces (collagen). The pieces are held together by suture. Pictures are below in **Figure 2.1** and **Figure 2.2**.

Figure 1.1 Angio-Seal™ VIP VCD Sealing the Arterial Puncture

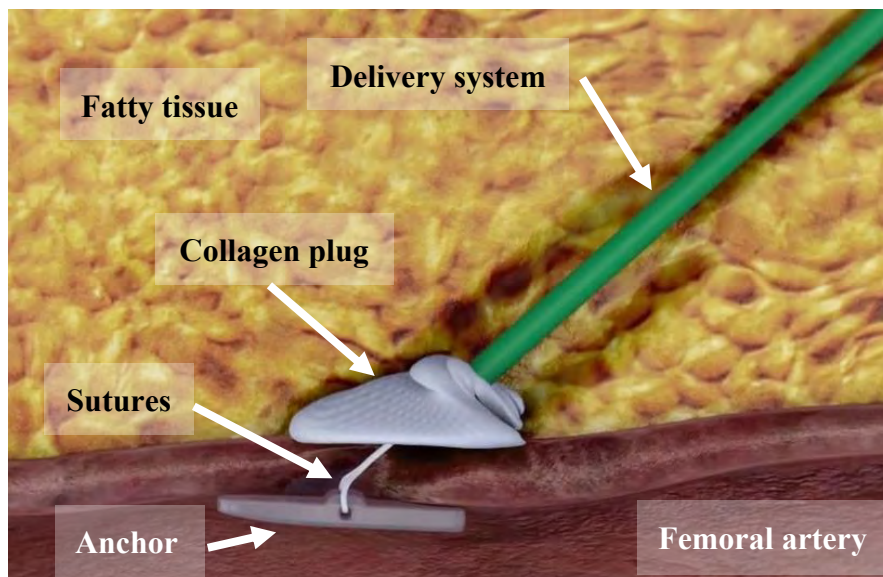
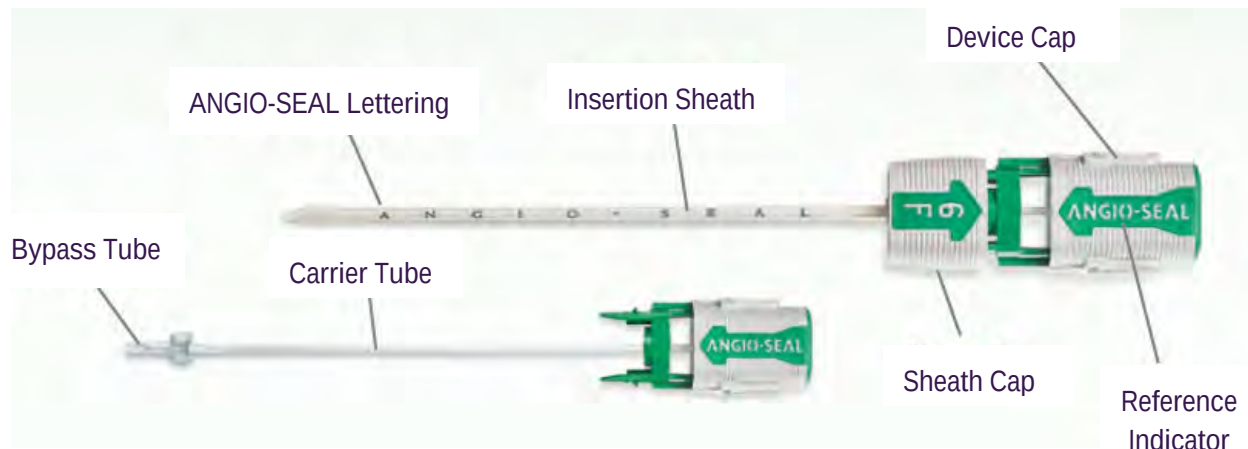
Figure 1.2 Angio-Seal™ VIP VCD

1.2.2 Raw Materials & Components

All the parts left in the body will dissolve over time.

Table 2.1: Angio-Seal™ VIP VCD Implantable and Non-Implantable Components

Unit	Part Name		Material
Bioabsorbable Unit *	Anchor	50:50 poly D,L lactide glycolide copolymer	
	Suture	Polyglycolic acid Polycaprolate / polyglycolic acid copolymer coating	
	Collagen Plug	Bovine hides	
Delivery System**	Carrier Tube	Nylon or High-Density Polyethylene (HDPE)	
	Tamper Tube	Polyurethane (6F) Polyethylene (8F)	
	Insertion Sheath	Insertion Sheath Tubing	Polyether Block Amine (PEBAX) MDX/Hexane Solution
		Molding	Acrylonitrile Butadiene Styrene (ABS)
		Sleeve	Polyurethane (8F only)
		Hemostasis Seal	Silicone/ Silicone Lubricant
		Insertion Sheath Cap	Acrylonitrile Butadiene Styrene (ABS)
Arteriotomy Locator	Arteriotomy Locator	High Density Polyethylene (HDPE)	

		Silicone lubricant
Accessories	Guidewire with J-straightener 6Fr: 0.035" OD, 70 cm length guidewire 8Fr: 0.038" OD, 70 cm length guidewire	Guidewire: 304 Stainless Steel J-Straightener: Polypropylene

* Implantable portion of device

**Angio-Seal™ VIP Vascular Closure Device Delivery System consists of the delivery device containing the bioabsorbable unit and the insertion sheath.

1.2.3 Accessories Intended to Be Used in Combination with the Device

Angio-Seal™ VIP VCD does not require additional accessories for the device to function as intended.

1.3 Risks and Warnings

Contact your healthcare professional if you believe that you are experiencing side effects related to the device or its use or if you are concerned about risks. This document is not intended to replace a consultation with your healthcare professional if needed. Your doctor has considered the risks and decided that this device is appropriate for your procedure.

1.3.1 Risks

Remaining Risks

Harms associated with the use of the Angio-Seal™ VIP VCD as a result of residual risks are listed in the Table below. The occurrence rate is based on data collected from January 2018 through April 2023.

Table 1.2: Angio-Seal™ VIP VCD Adverse Events/Complications

Potential Harms	N = number of adverse events	Occurrence Rate (%) N/total sales	Mitigating Factors
Allergic Reaction	4	0.0000	Actions and plans to reduce risk as low as possible have been addressed by product design features, manufacturing guidelines, product labeling, and physician training.
Aneurysm	2	0.0000	
AV Fistula	1	0.0000	
Blood Loss/Bleeding	1774	0.0216	
Ecchymosis	4	0.0000	
Embolism	23	0.0003	
Foreign Body Reaction	54	0.0007	
Hematoma	132	0.0016	
Hemorrhage	13	0.0002	
Infection	11	0.0001	

Inflammatory Reaction	9	0.0001
Pain	15	0.0002
Patient discomfort	640	0.0078
procedure delay	1165	0.0142
Pseudoaneurysm	45	0.0005
Retroperitoneal Bleed	22	0.0003
Thromboembolism	14	0.0002
Thrombosis	18	0.0002
Vessel Occlusion	136	0.0017
Vessel Perforation	3	0.0000
Vessel Tissue Dissection/Laceration	8	0.0001

1.3.2 Warnings and Precautions

- Follow doctor's orders.
- Call your doctor if you have bleeding, pain, or signs of infection. Signs of infection are redness, fever, discharge an/or warmth at the site. Signs of blood loss (hematoma) are swelling and lump in groin.
- Carry the Patient Information Card for the next 90 days.

1.3.3 Field Safety Corrective Actions

There were no Field actions and recalls for Angio-Seal™ VIP VCD between January 1, 2018, and April 30, 2023.

1.4 Clinical Evaluation and Post-Market Clinical Follow-Up

1.4.1 Summary of Data from Animal Testing

Several animal studies were conducted to support FDA approval and CE mark certification of Angio-Seal™ VCD. Pre-clinical evidence is based on equivalence to other versions of the Angio-Seal™ Device. Details about these studies are provided in the table below:

Table 2.3: Summary of Data from Animal Testing

Report Number	Device	Time to Hemostasis (minutes)	Hematoma (%)	Deployment Force (lbf)
VAL-90017642	Angio-Seal™ Evolution VCD 6F (6 punctures)	0.83 ± 2.04	17%	N/A
VAL-90017642	Angio-Seal™ Evolution VCD 8F (6 punctures)	0.67 ± 1.63	33%	N/A
VAL-90017642	Angio-Seal™ VIP VCD 6F (3 punctures)	0 ± 0	0%	N/A
VAL-90017642	Angio-Seal™ VIP VCD 8F (3 punctures)	0.83 ± 1.44	0%	N/A
VAL-90018838	Angio-Seal™ Evolution VCD 8F (4 punctures)	1.25 ± 2.50	25%	3.94 ± 0.8
VAL-90018838	Angio-Seal™ VIP VCD 8F (4 punctures)	2.88 ± 4.52	0%	2.98 ± 0.5

1.4.2 Summary of Data from Pre-Market and Post-Market Clinical Investigations

The table below includes a list of clinical and post-market studies. Clinical evidence is not based on equivalence to other versions of the Angio-Seal™ Device.

Table 2.4: Summary of Clinical Studies

First in Human Pre-market Angio-Seal – original iteration (HPCD) Exploratory Clinical Study (8 French)	Iteration of Angio-Seal	MDCG 2020-06 Rank
Kensey Nash Hemostatic Puncture Closure Device (HPCD™) Phase I (G910028) (1992) 1. US Phase 1 2. US Phase 1 Extended 3. European Phase 1 4. European Phase 1 Extended	Angio-Seal Hemostatic Closure Devices (original iteration)	Rank 2
First in Human Pre-market Angio-Seal – original iteration (HPCD) Exploratory Clinical Study (6 French)	Iteration of Angio-Seal	MDCG 2020-06 Rank
6 French Pilot Study	Angio-Seal Hemostatic Closure Devices (original iteration)	Rank 4
Pre-market confirmatory Studies - original iteration (HPCD)	Iteration of Angio-Seal	MDCG 2020-06 Rank
US Phase II European Phase II Early Ambulation in Diagnostic Angiography Patients	Angio-Seal Hemostatic Closure Devices (original iteration)	Rank 4
Pre-Market Approval Studies (Leading to Angio-Seal VIP VCD)	Iteration of Angio-Seal	MDCG 2020-06 Rank
Clinical Evaluation of the Angio-Seal Self-Tightening Slipknot (STS) Hemostatic Puncture Closure Device (IDE G990306)	Angio-Seal STS	Rank 2
Angio-Seal STS Early Ambulation and Discharge Study	Angio-Seal STS	Rank 2
Completed Post-Market Clinical Studies	Iteration of Angio-Seal	MDCG 2020-06 Rank

Clinical Study of Angio-Seal™ VIP Vascular Closure Device following Diagnostic and/or Interventional Endovascular Procedures (CP0901) (2009)	Angio-Seal™ VIP VCD	Rank 2
SEAL PM Post-Market Study of the St. Jude Medical Angio-Seal™ VIP Vascular Closure Device CV-13-008-EU-VC (2013)	Angio-Seal™ VIP VCD	Rank 2
Completed Post Market Clinical Follow Up Activities	Iteration of Angio-Seal	MDCG 2020-06 Rank
20210293 – Angio-Seal VIP Vascular Closure Device Clinical Survey Report (2020)	Angio-Seal™ VIP VCD	Rank 4
20210298 – Guidewire Clinical Survey Report (2020)	Angio-Seal™ VIP VCD	Rank 4
Post-Market Animal Testing	Iteration of Angio-Seal	MDCG 2020-06 Rank
Angio-Seal® Vascular Closure Device Biological Evaluation & Testing Report 17-3-121 (2018)	Angio-Seal™ VIP VCD Angio-Seal™ Evolution	Rank 12
Animal Study Report, Evolution CER044, Chronic Animal Study 93d, APS #90014675 (2008)	Angio-Seal™ VIP VCD Angio-Seal™ Evolution	Rank 11
Animal Study Report, Evolution CER056, Acute Animal Study at APS #90017642 (2008)	Angio-Seal™ VIP VCD Angio-Seal™ Evolution	Rank 11
Animal Study Report, Evolution CER056, Acute Animal Study at APS #90018838 (2008)	Angio-Seal™ VIP VCD Angio-Seal™ Evolution	Rank 11
Planned / Ongoing Post-Market Clinical Studies	Iteration of Angio-Seal	MDCG 2020-06 Rank
The Angio-Seal Close: Post-market Clinical Investigation of the Angio-Seal™ VIP Vascular Closure Device (VCD): A Prospective, Multi-Center, Observational Study in Europe (T139E4)	Angio-Seal™ VIP VCD	Rank 3 (ongoing)

The Angio-Seal™ VIP belongs to the Angio-Seal™ Family of devices which have been in Europe since 1994 and in the U.S. since 1996. The objective of the clinical studies was to confirm the safety and effectiveness of Angio-Seal™.

Table 2.5: First in Human and Pre-market Studies

First in Human Pre-market Angio-Seal – original iteration (HPCD) Exploratory Clinical Study (8 French)					
Exploratory Studies - US Phase 1 - US Phase 1 Extended - European Phase 1 - European Phase 1 Extended					
Name of Study / Year: US Phase I Kensley and Nash Hemostatic Puncture Closure Device – Phase I Clinical Study 1992 Device Iteration: Angio-Seal™ Hemostatic Puncture Closure	Study Design Type: Multi-centered non-randomized exploratory study	Country: United States Follow Up time points: - 24-hours post op - 7-day post op - 30-day post op - 6-week post op - 60-day post op	Study Endpoints: deployment success, device failure, immediate hemostasis, complications, collagen absorption via US Performance Outcomes Successful device deployment 19/20	Potential Gaps: 6 Fr devices are not used (both for diagnostic and interventional procedure) 8 Fr devices are not used for interventional procedure Time to hemostasis definition “This is defined	Justification of Gaps/ How gaps were/are being addressed: Due to the stage of development, an exploratory study is appropriate. In order to best assess device efficiency and general

First in Human Pre-market Angio-Seal – original iteration (HPCD) Exploratory Clinical Study (8 French)					
Device (HPCD) – original iteration			Device failures 1/20 -1 device did not anchor in the vessel and simply pulled out with the introducer; hemostasis was achieved via manual pressure in this patient. Patient was removed from follow up Immediate hemostasis 11/19 - 5 devices produced hemostasis in under 3 minutes - 1 device produced hemostasis in 5 minutes with gentle manual pressure - 1 device produced hemostasis in 10 minutes with accompanying manual pressure - 1 device produced hemostasis with accompanying manual pressure after prolonged slow seeping - 1 device did not anchor in the vessel and simply pulled out with the introducer; hemostasis was achieved via manual pressure in this patient. Patient was removed from follow up. (Same device listed under device failure)	as the time elapsed from device deployment (cutting of the suture) to complete cessation of arterial bleeding” exists but does not align with current definition Time to ambulation not defined or collected Exclusion criteria beyond the one defined in IFU as special patient population -Exclusion criteria excludes patients commonly known to have endovascular procedures that would utilize femoral artery closure via VCD. -Exclusion criteria excluded patients with any vascular compromise of the leg to be involved including previous amputation, peripheral vascular surgery, interventional procedure, clinically significant peripheral vascular disease, or presence of a femoral artery bruit. Guidewire included as a device accessory was not evaluated as part of the study.	safety, the variables and risk factors in the subject population were limited to avoid biases and skew the data. Excluding subjects who were obese, had a compromised femoral artery(ies), or evidence of peripheral artery disease removes any factors that may be misinterpreted as a device failure rather than an anatomy not suited for use of the device. To address the time to hemostasis for the 6Fr diagnostic and interventional procedures, Clinical Study of Angio-Seal™ VIP Vascular Closure Device following Diagnostic and/or Interventional Endovascular Procedures (CP0901) (NCT01393041) and SEAL PM Post-Market Study of the St. Jude Medical Angio-Seal™ VIP Vascular Closure Device CV-13-008-EU-VC (2013) NCT01858636 were performed on the Angio-Seal VIP subject device. For the definition of time to hemostasis: the sponsor has created a standardized
Published/un-published: Un-published	Inclusion Criteria - Over 18 years old - Able to give written informed consent - Clinically indicated for invasive procedure using Daig 8 French introducer sheath in the femoral artery.				
MDCG 2020-06 Rank: Rank 2	Exclusion Criteria - Known allergy to the materials in the device. - Obesity which would result in an unusual approach to the femoral artery.	Patient Demographics: -Age range 41-71 -Mean age 59 -Males 16 -Females 4	Safety / Adverse Event Outcomes - Hematoma 2 - 1 measured 3 X 5 cm - 1 measured 15 cm diameter		

First in Human Pre-market Angio-Seal – original iteration (HPCD) Exploratory Clinical Study (8 French)					
Study Type Per ISO 14155: Exploratory (initial feasibility and acute safety of the device)	- Severe acute systemic disease or terminal illness. - Any infection. - Bleeding diathesis. - Any vascular compromise of the leg to be involved including previous amputation, peripheral vascular surgery, interventional procedure, clinically significant peripheral vascular disease, or presence of a femoral artery bruit. - Pregnancy - Patients on inotropic or intra-aortic balloon pump support. - Emergency catheterization procedures. - Use of thrombolytic agents within 24 hours prior to, during, or after the catheterization procedure. - Systolic blood pressure of greater than 180 in either arm. - Ankle/brachial index less than 0.9 in either leg.	Procedural characteristics - Diagnostic angiography 17 - PTCA 3	- Both resolved spontaneously - Ecchymosis 8 - Infections 0 - Other complications 0		definitions of time to hemostasis, time to ambulation and oozing. These definitions will be utilized going forward and are included in the ongoing PMCF Angio-Seal™ Close Registry. This registry also includes data collection on time to ambulation and the guidewire that is included with the Angio-Seal™ VIP VCD as a device accessory. The guidewire included as a device accessory was also evaluated in a Rank 4 High Quality Clinical Survey: 0.035" and 0.038" Guidewire Clinical Survey Report.
Study Objective To examine the initial feasibility and acute safety of the device (HPCD) in 20 patients for whom unilateral femoral access was indicated for therapeutic or diagnostic reasons			Imaging results (US) - 24-hours post op: anchors well placed and flush against intra-arterial wall - 7-day post op: anchors less echogenic but still detectable - 30-day post op: (19 patients) anchor visible in only 8 patients, not detectable in 11 patients - 6-week post op: (1 patient) evidence of anchor, but no additional follow up on this patient. - 60-day post op: (15 patients) minimal evidence of anchor in 2 patients		
Name of Study / Year: US Phase I Extended Kensey Nash Hemostatic Puncture Closure Device IDE No. G910028 Investigational Plan for Extended Phase I study 1992 Device Iteration: Angio-Seal™ Hemostatic Puncture Closure Device (HPCD) – original iteration	Study Design Type: Multi-centered non-randomized exploratory study Inclusion criteria - Over 18 years - Under 80 years - Able to provide written informed consent - Clinically indicated for invasive procedure involving femoral arterial access	Country: United States Follow Up time points: Changes that were requested and approved included - Elimination of 7 day and 1 month US studies in favor of 1 US post-HPCD placement, with 90-day post placement US available at investigators discretion - Rationale was accepted	Study Endpoints: deployment success, device failure, immediate hemostasis, complications, collagen absorption via US Performance Outcomes: Successful anchor and device deployment 9/10 (90%) Hemostasis - Immediate hemostasis achieved 8/9 (88%) - Hemostasis within 1-minute 1/9(11%)	Potential Gaps: - US Phase I Extended - Study was voluntarily suspended by the sponsor after enrollment of 10 patients due to failure in of an anchor in a study in USA	Justification of Gaps/ How gaps were/are being addressed: Per clinical study report, anchor detachment was traced back to certain aspects of the manufacturing process that could have introduced weakness into the anchor. This was addressed by eliminating the aspects of the manufacturing processes by a

First in Human Pre-market Angio-Seal – original iteration (HPCD) Exploratory Clinical Study (8 French)					
Published/un-published: Un-published			- Seeping upon removal of tamper – 1/9 (11%) (Resolved with pressure applied with an IV bag and dressing) - Seeping for 30 & 60 minutes post-Angio-Seal placement – 2/9 (22%) (Resolved with gentle pressure) - Anchor performance -1/10 patients, 2 devices attempted and did not anchor in the vessel and simply pulled out the introducer; 3rd device was introduced and deployed; hemostasis produced; required digital pressure for 2 minutes		change in the anchor design. Due to study suspension, the interpretation of the data will not be done.
MDCG 2020-06 Rank: Rank 2	Exclusion Criteria: - Allergy to any of the materials used in the device. - Any severe acute systemic disease or terminal illness. - Any evidence of infection, including groin infection. - Any bleeding diathesis (except use of aspirin or heparin, both of which are allowed). - Clinically significant peripheral vascular disease, previous percutaneous procedure (diagnostic or therapeutic, using the leg involved) in the 7 days prior to use of the HPCD, or previous percutaneous procedure (diagnostic or therapeutic, using the leg involved) earlier than 7 days prior to the HPCD where there is evidence of remaining	Patient Demographics: - Age range – 54-78 years - Mean age 65 years - Males 7 - Females 3 Procedural characteristics: -Diagnostic angiography 9 -PTCA 1	Safety / Adverse Event Outcomes: - Hematoma 1 -1, 3 cm hematoma and bruising, treated with manual pressure - Ecchymosis 1 - Infections 0 - Other complications 0 - Device failures 0 -1 patient lost to follow up due to death, 3 months post-placement, which was due to acute myelogenous leukemia, unrelated to Angio-Seal device Imaging: - US completed within 24 hours post-placement -All anchors well-placed and flushed		

First in Human Pre-market Angio-Seal – original iteration (HPCD) Exploratory Clinical Study (8 French)					
	hematoma, ecchymosis, palpable mass or tenderness. - Pregnancy. - Patients on inotropic or balloon pump support. Emergency catheterization procedures. - Use of thrombolytic agents within 24 hours prior to or during the catheterization procedure. Systolic blood pressure >180 mm Hg in either arm.		against the arterial wall 60 day follow up -5 patients detected only minimal evidence of the anchor in one patient 9 weeks -1 patient normal, showing no evidence of anchor 5 months - 1 patient, normal US with no anchor evident US studies revealed no evidence of stenosis, thrombosis, thrombo-embolism, changes in flow patterns, or changes in arterial architecture.		
Study Type Per ISO 14155: Exploratory (initial feasibility and acute safety of the device)					
Study Objective The device being tested is the Kensey Nash Hemostatic Puncture Closure Device (abbreviated HPCD™). This is a continuation of the Phase I clinical study, termed “Extended Phase I,” will assess the efficacy and safety of the Hemostatic Puncture Closure Device (HPCD) in humans.					
Name of Study / Year: European Phase I European Phase I Extended (Data reported combined)	Study Design Type: Multi-centered non-randomized exploratory study	Country: Netherlands Follow Up time points: - 1 week - 1 month	Study endpoints: deployment success, device failure, immediate hemostasis, complications, collagen absorption via US	Potential Gaps 6 Fr devices are not used (both for diagnostic nor interventional procedure) 8 Fr devices are not used for interventional procedure	Justification of Gaps/ How gaps were/are being addressed: To address the time to hemostasis for the 6Fr diagnostic and interventional

First in Human Pre-market Angio-Seal – original iteration (HPCD) Exploratory Clinical Study (8 French)					
Device Iteration: Angio-Seal™ Hemo-static Puncture Closure Device (HPCD) – original iteration Published/un-published: Un-published		- 3 months	Performance Outcomes Primary hemostasis 39/55 (70.09%) Hemostasis < 2 minutes 6 Hemostasis 2-5 minutes 5 Hemostasis 10 minutes 1 Blood seeped (slow trickle of blood) > 12 hours 1 Patient did not require transfusion Required manual pressure to be applied 2 Non anchoring of device or device failure 3	Time to hemostasis definition exists but does not align with current definition Time to ambulation not defined or collected There were no exclusion criteria beyond the one defined in IFU as special patient population Guidewire included as a device accessory was not evaluated as part of the study.	procedures, Clinical Study of Angio-Seal™ VIP Vascular Closure Device following Diagnostic and/or Interventional Endovascular Procedures (CP0901) (NCT01393041) and SEAL PM Post-Market Study of the St. Jude Medical Angio-Seal™ VIP Vascular Closure Device CV-13-008-EU-VC (2013) NCT01858636 were performed on the Angio-Seal VIP subject device. The definition of time to hemostasis, the sponsor has created a standardized definitions of time to hemostasis, time to ambulation and oozing. These definitions will be utilized going forward and are included in the ongoing PMCF Angio-Seal™ Close Registry. This registry also includes data collection on time to ambulation and the guidewire that is included with the Angio-Seal™ VIP VCD as a device accessory. The guidewire included as a device accessory was also evaluated in the Clinical Survey *0.35 and 0.38
MDCG 2020-06 Rank: Rank 2 Study Type Per ISO 14155: Exploratory (initial feasibility and acute safety of the device) Study objective: To examine the initial feasibility and acute safety of the device (HPCD) in a small group of patients for whom unilateral femoral access was indicated for therapeutic or diagnostic reasons.	Inclusion Criteria -Over 18 years old -Able to give informed consent -Clinically indicated for an invasive procedure using an 8F introducer sheath in femoral artery Exclusion Criteria: - Allergy to the materials in the device. - Obesity which would result in an unusual approach to the femoral artery (Phase I only). -Severe acute systemic disease or terminal illness. - Any infection. - Bleeding diathesis. - Any vascular compromise of the leg to be involved including previous amputation, peripheral vascular surgery, interventional procedure, clinically significant peripheral vascular disease, or presence of a	Patient Demographics: n= 55 patients Male 37 Female 18 Male age range 42-74 Male mean age 56.8 Female age range 45-72 Mean female age 61.1 aPTT > 300 seconds 46 aPTT 200-300 seconds 1 aPTT 100-200 seconds 1 aPTT < 100 seconds 6 Procedural characteristics Diagnostic angiography 20 Percutaneous transluminal coronary angioplasty (PTCA) 35	Safety / Adverse Event Outcomes Secondary bleeding within 30 minutes 8 -2 required pressure to be applied to produce hemostasis Secondary bleeding stopped > 30 minutes 19 -8 required pressure to produce hemostasis Ecchymosis 24 Palpable mass at puncture site 7-day post-placement 9 -All resolved in 30 days Hematomas 2 Ambulated 4 hours post-placement 26 Ambulated > 7 hours 9 1 ambulation unknown 3 patients ambulation time N/A		

First in Human Pre-market Angio-Seal – original iteration (HPCD) Exploratory Clinical Study (8 French)					
	femoral artery bruit (Phase I) or clinically significant peripheral vascular disease, previous percutaneous procedure using the leg involved in the 7 days prior to use of the Angio-Seal, or previous percutaneous procedure earlier than 7 days prior to the Angio-Seal where there is evidence of remaining hematoma, ecchymosis, palpable mass or tenderness (Phase I Extended) - Pregnancy - Patients on inotropic or intra-aortic balloon pump support - Emergency catheterization procedures. - Use of thrombolytic agents within 24 hours prior to, during, or after the catheterization procedure. - Systolic blood pressure of greater than 180 in either arm. - Ankle/brachial index less than 0.9 in either leg.		Device failed to anchor in the vessel 2 -Manual pressure was applied to produce hemostasis Positioning suture broke 2 -1 related to operator error Collagen protrusion from skin after device placement 1 Imaging All US studies unremarkable No cases of intraluminal thrombosis, residual stenosis at Angio-Seal site Change in flow spectrum, spectral broadening noted in some patients indicating slightly more turbulent flow, resolved by 1 week follow up in majority of patients; investigator did not consider this a significant finding.		Guidewire Survey Report”

Table 2.6: First in Human Pre-market Studies

First in Human Pre-market Angio-Seal – original iteration (HPCD) Exploratory Clinical Study (8 French)					
Name of Study / Year:	Study Design Type:	Country:	Study Endpoints:	Potential Gaps:	Justification of Gaps/ How gaps were/are being addressed:
French Pilot Study 1996	Single center pilot nonrandomized controlled clinical study	France	Time to hemostasis, time to ambulation, time to discharge, and complications (device related or otherwise)	8 Fr devices are not used 6 Fr device for interventional procedure was not used Time to hemostasis definition exists but does not align with current definition	The 6French Angio-Seal device was not marketed at the time of this study. To address the time to hemostasis for the 6Fr diagnostic and
Device Iteration: 6Fr Angio-Seal device, Sherwood Medical Co (dba The Kendall Company, manufactured by Kensey Nash)	Historical manual pressure group utilized for comparisons	Follow Up time points: 7 to 14 days	Performance Outcomes Successful deployment 47/50 (94%) - Procedure was aborted or access was lost due to difficulty inserting the sheath and arteriotomy locator during sheath exchange or		

			difficulty in advancing the device through the sheath in 2 patients - 1 device non-deployment due to physician inserting the device into the sheath upside down. Time to hemostasis - Immediate (time to hemostasis =0) 98% - 30 minutes 87% Time to discharge Device deployment success in training patients 13/16 (81%) -Prior to patient enrollment in study 16 patients were training patients.	Time to ambulation does not align with current definition There were no exclusion criteria beyond the one defined in IFU as special patient population Guidewire included as a device accessory was not evaluated as part of the study.	interventional procedures, Clinical Study of Angio-Seal™ VIP Vascular Closure Device following Diagnostic and/or Interventional Endovascular Procedures (CP0901) (NCT01393041) and SEAL PM Post-Market Study of the St. Jude Medical Angio-Seal™ VIP Vascular Closure Device CV-13-008-EU-VC (2013) NCT01858636 were performed on the Angio-Seal VIP subject device. The definition of time to hemostasis, the sponsor has created a standardized definitions of time to hemostasis, time to ambulation and oozing. These definitions will be utilized going forward and are included in the ongoing PMCF Angio-Seal™ Close Registry. This registry also includes data collection on time to ambulation and the guidewire that is included with the Angio-Seal™ VIP VCD as a device accessory. The guidewire included as a device accessory was also evaluated in the Clinical Survey “2021 Clinical Survey”
<u>MDCG 2020-06 Rank:</u> Rank 4 <u>Published/un-published:</u> Un-published <u>Study Type Per ISO 14155:</u> Exploratory (initial feasibility and acute safety of the device) <u>Study Objective</u> Pilot study to explore the feasibility and acute safety of the device (HPCD) in a small group of patients for whom unilateral femoral access using 6Fr device was indicated diagnostic reasons.	<u>Inclusion criteria:</u> - Over 18 years of age. - Able and willing to give informed consent. - Clinically indicated for diagnostic angiography involving access through the femoral artery. - Patient expected to be physically able to ambulate for 100 feet or five minutes following discharge from the cath lab, completion of the hemostasis procedure, and appropriate bedrest as determined by the protocol and expected to be available for follow-up at 7-14 days as assessed by the investigator.	<u>Patient Demographics:</u> N= 50 patients <u>Procedural characteristics</u> -Diagnostic angiography 50	<u>Safety / Adverse Event Outcomes</u> Including Device Malfunction: Any Complication, AS Group: 14/200 (7.0%); MP Group: 2/93 (2.2%); No Complication, AS: 186 (93.0%); MP: 91 (97.8%); Excluding Device Malfunction, AS: 12 (6.0%); MP: 2 (2.2%); No Complication, AS: 188 (94.0%); MP: 91 (97.8%) <u>Complications (AS vs MP):</u> Vascular repair, 1 (0.5%) vs 0 (0%), P=1.000; Bleeding requiring transfusion, 0 (0%) vs 0 (0%) Infection + hospital, 0 (0%), both groups Device malfunction, 2 (1.0%) in AS group Hematoma (≥6 cm), 1 (0.5%) vs 0 (0%), P=1.000		

	- Females are eligible for the study provided that the investigator has determined they are not pregnant or lactating, are surgically sterile, post-menopausal, or have a negative pregnancy test. Exclusion Criteria: - Known allergy to any of the materials used in the device - Evidence of systemic bacterial or cutaneous infection, including groin infection - Potential coagulopathy as indicated by platelet count < 75,000 - Use of systemic thrombolytic agents within 24 hours prior to or during the catheterization procedure which cause the concentration of fibrinogen to be < 100 mg/dl. - Sheath used for invasive procedure in place for greater than 36 hours - Sheath size used for invasive procedure larger than 6 French - Use of any other unapproved investigational drug or device in the last 30 days - Suspected double wall puncture at the site to be closed with the Angio-Seal device. - If a palpable hematoma (> 2 cm in diameter) is observed at the end of the procedure, the patient should be excluded from the study.		Hematoma (<6 cm), 3 (1.5%) vs 0 (0%), P= 0.554 Late bleeding 7 (3.5%) vs 2 (2.2%) P=0.724 CI: 1.6 [0.4-7.7] Pseudoaneurysm 1 (0.5%) vs 0 (0%), P= 1.000 AV Fistula 1 (0.5%) vs 0 (0%), P=1.000 DVT 1 (0.5%) vs 0 (0%), P= 1.000 No major complication, 199 (99.5%) vs 93 (100%), P= 1.000		

	- If there is any indication that the puncture has been made in the profunda femoris or the bifurcation of the common femoral artery, the patient should be excluded from the study.				

Table 2.7: Pre-market Confirmatory Studies

Pre-market confirmatory Studies - original iteration (HPCD)					
Confirmatory Studies 1. US Phase II (Reference: PMA P930038 Volume 4, Page VI-32) 2. European Phase II (Reference: PMA P930038 Volume 4, page VI-80) 3. Early Ambulation in Diagnostic Angiography Patients (Reference: Supplemental Change: PMA P930038/S029)			30 September 1996 –Angio-Seal HPCD FDA conditional approval, Approval was conditional upon a post approval study to be conducted to examine the incidence of infection and vascular compromise associated with Angio-Seal device use. Follow up US information should be analyzed by a blinded core laboratory for objective interpretation of data.		
Name of Study /	Study Design Type:	Country:	Study Endpoints:	Potential Gaps:	Justification of Gaps/ How gaps were/are being addressed:
Year: US Phase II	Prospective, randomized controlled study	United States	Safety Endpoints: Hemostasis Time (rapid or prolonged); Late Bleeding or Oozing; Hematoma, Thrombosis, false aneurysm or arteriovenous fistula; Infection at puncture site; Allergic reaction to the device; and Device failure	8 Fr devices not used for interventional procedures. 6 Fr devices are not used (not marketed until 1999) Definition of Time to Hemostasis is included but is not aligned with our current definition. Time to Ambulation was not defined or collected. Exclusion criteria:	The 6 French Angio-Seal device was not marketed at the time of this study. To address the time to hemostasis for the 6Fr diagnostic and interventional procedures, Clinical Study of Angio-Seal™ VIP Vascular Closure Device following Diagnostic and/or Interventional Endovascular Procedures (CP0901) (NCT01393041) and SEAL PM Post-Market Study of the St. Jude Medical Angio-Seal™ VIP Vascular Closure Device CV-13-008-EU-VC (2013) NCT01858636 were performed on the Angio-Seal VIP subject device. To address the definition of time to hemostasis, the
Device Iteration: Angio-Seal™ Hemostatic Puncture Closure Device (HPCD) – original iteration			Follow Up time points: 2 months Performance End-points: Time to hemostasis (the time required for the bleeding from the puncture to stop), compression time (the time manual pressure was applied to produce hemostasis) Performance Outcomes: - Time to Hemostasis: Angio-Seal group: 2.5±9.2 mins on average		

			- Control group: 16.2±12.6 mins (p<0.0001) <u>Compression Time:</u> - Angio-Seal group: 2.5±9.2 mins on average - Control group: 20±12.6 mins (p<0.0001) - 5-mins post-sheath removal, 86.1% Angio-Seal group patients while only 6.5% control group patients achieved hemostasis	Patients older than 80 were excluded Patients excluded if sheath placed for >36 hours, which is not on current list of Special Populations. Guidewire included as a device accessory was not evaluated as part of the study.	sponsor has created a standardized definitions of time to hemostasis, time to ambulation and oozing. These definitions will be utilized going forward and are included in the ongoing PMCF Angio-Seal™ Close Registry. This registry also includes data collection on time to ambulation and the guidewire that is included with the Angio-Seal™ VIP VCD as a device accessory. The guidewire included as a device accessory was also evaluated in the Clinical Survey
<u>MDCG 2020-06 Rank:</u> Rank 2	<u>Inclusion Criteria:</u> - Male and female subjects over age 18 and under age 80 - Patients were required to give written informed consent	<u>Patient Demographics:</u> N= 311	<u>Safety / Adverse Event Outcomes</u> Complication Rates - Angio-Seal Group: 17.8% - Control Group: 22.2 (p=0.334) Non-bleeding Complication Rates: - Angio-Seal Group: 8.9% - Control Group: 19.6 (p=0.011) The most frequent complication was bleeding in the Angio-Seal group while, it was long compression time (≥30 mins) in the control group. However, the severity of complications was similar in both the groups.		
<u>Study Type Per ISO 14155:</u> Confirmatory	- Clinically indicated for an invasive procedure involving access through an 8F or smaller puncture in the femoral artery	<u>Procedural characteristics:</u> Angio-Seal (Device Group) 157 patients. Manual Pressure (Control Group) 154 patients			
<u>Study objective</u> Phase II clinical trials were designed to study the safety and efficacy of the Angio-Seal in a larger number of patients in up to 12 investigational sites. Phase II clinical trials compared the Angio-Seal to the current standard of care in arterial puncture hemostasis which is the application of manual pressure to the puncture site. Efficacy was	<u>Exclusion Criteria:</u> - Known allergy to any of the materials used in the device - Severe acute non-cardiac systemic disease or terminal illness - Evidence of systemic bacterial or cutaneous infection, including groin infection - Possible bleeding diathesis as indicated by platelet count <125,000; patients on anti-coagulant therapy were excluded if prothrombin time (PT) >15 seconds				

evaluated in two categories: time to hemostasis and compression time. The safety of the device was evaluated by evaluating any individual complications and the overall complication rate.	(heparin and ASA allowed). - Clinically significant peripheral vascular disease of the proposed artery as indicated by the following: a. Foot pulses absent in both feet b. Femoral bruit c. Claudication d. Prior vascular surgery on the femoral artery e. Ankle/Brachial index less than 0.9 - Previous arterial access (for diagnostic or therapeutic purposes, using the proposed artery) in the 7 days prior if duplex ultrasound of the site is abnormal - Women of child-bearing potential - Inotropic or balloon pump support - Use of thrombolytic agents within 24 hrs prior to or during the catheterization procedure or history of recent (within 7 days) thrombolysis and fibrinogen <100 mg/d/dL on day of the study - Baseline systolic blood pressure greater than 180 mm Hg in either arm - Sheath used for invasive procedure in place for longer than 36 hrs				
Name of Study / Year: European Phase II	Study Design Type: Prospective randomized controlled study	Country: Europe	Study Endpoints: Time to hemostasis, incidence of any groin-site complications such as bleeding, hematoma, ecchymosis or other abnormalities; and visual evidence of collagen plug and anchor via US.	Potential Gaps: 6 Fr devices are not used (not marketed until 1999)	Justification of Gaps/ How gaps were/are being addressed: To address the time to hemostasis for the 6Fr diagnostic and interventional procedures, Clinical Study of Angio-Seal™ VIP Vascular Closure Device following Diagnostic and/or Interventional Endovascular Procedures (CP0901)
Device Iteration: Angio-Seal™ Hemostatic Puncture Closure Device (HPCD)		Follow Up time points: - 2 months post placement - 90 days at the discretion of the investigator		Definition of Time to Hemostasis is included but is not aligned with our current definition.	

– original iteration			Performance Outcomes: - A shorter average time to hemostasis in Angio-Seal group (1.13 min) when compared to the manual pressure group (7.51 min) - Ambulation: 42/55 Angio-Seal group patients (76%) ambulated within four hours while in manual pressure group, only 1/52 patients ambulated within six hours. Only one minor complication was associated with early (<6 hrs) ambulation	Time to Ambulation was collected, but not clearly defined. Exclusion Criteria: Patients older than 80 were excluded Patients excluded if sheath placed for >36 hours, which is not on current list of Special Populations.	(NCT01393041) and SEAL PM Post-Market Study of the St. Jude Medical Angio-Seal™ VIP Vascular Closure Device CV-13-008-EU-VC (2013) NCT01858636 were performed on the Angio-Seal VIP subject device. To address the definition of time to hemostasis, the sponsor has created a standardized definitions of time to hemostasis, time to ambulation and oozing. These definitions will be utilized going forward and are included in the ongoing PMCF Angio-Seal™ Close Registry. This registry also includes the guidewire that is included with the Angio-Seal™ VIP VCD as a device accessory. The guidewire included as a device accessory was also evaluated in the Clinical Survey
MDCG 2020-06 Rank: Rank 2 Study Type Per ISO 14155: Confirmatory Study objective The objective of the study is to assess the efficacy and safety of the Angio-Seal compared to the current standard of care in puncture hemostasis which is the application of manual compression to the arterial puncture site. The safety of the device was evaluated by examining the performance of the device versus the performance of manual compression and by evaluating any individual complications and overall complications.	Inclusion Criteria: - Male and female subjects over age 18 and under age 80 - Patients were required to give written informed consent - Clinically indicated for an invasive procedure involving access through an 8F or smaller puncture in the femoral artery Exclusion Criteria: - Known allergy to any of the materials used in the device - Severe acute non-cardiac systemic disease or terminal illness - Evidence of systemic bacterial or cutaneous infection, including groin infection - Possible bleeding diathesis as indicated by platelet count <125,000; patients on anti-coagulant therapy were excluded if prothrombin time (PT) >15 seconds (heparin and ASA allowed). - Clinically significant peripheral vascular disease of the proposed	Patient Demographics: N=107 Procedural Characteristics: Angio-Seal (Device Group) 55 patients; Manual Pressure (Control Group) 52 patients Procedure Details: Angio-seal Group (n=55): Diagnostic procedure, 33; therapeutic procedure, 22; Manual Pressure Group (n=52): Diagnostic procedure, 28; therapeutic procedure, 24;	Safety / Adverse Event Outcomes: - Complications reported in 14 patients/107 total patients: 12 in Angio-Seal group (12/55, 21.8%); 2 in manual pressure (2/52, 3.8%) group - The authors stated that the reported complications are all inclusive and are not broken down into adverse events or device failures and replacements. - All complications in both the groups were predominantly mild and none required intervention, surgical repair or blood transfusion. - The case in which the tamper was left in place for 7 days, II-MAAS-78, was clearly a failure to follow instructions for use. Along with, the two patients who had some bleeding when the device was deployed, II-MAAS-89 and II-MAAS-103, might not be considered complications since bleeding on deployment is expected in some patients; the bleeding, as time to hemostasis, is	Guidewire included as a device accessory was not evaluated as part of the study.	

	artery as indicated by the following: a. Foot pulses absent in both feet b. Femoral bruit c. Claudication d. Prior vascular surgery on the femoral artery e. Ankle/Brachial index less than 0.9 - Previous arterial access (for diagnostic or therapeutic purposes, using the proposed artery) in the 7 days prior if duplex ultrasound of the site is abnormal - Women of child-bearing potential - Ionotrophic or balloon pump support - Use of thrombolytic agents within 24 hrs prior to or during the catheterization procedure or history of recent (within 7 days) thrombolysis and fibrinogen <100 mg/dL on day of the study - Baseline systolic blood pressure greater than 180 mm Hg in either arm - Sheath used for invasive procedure in place for longer than 36 hrs		accounted for the efficacy analysis.		
Name of Study / Year: Early Ambulation in Diagnostic Angiography Patients STS Study in response to PMA July 5, 2001 – PMA P930038/S029 FDA approval letter <i>Changes occurred to the design of the device which triggered the need for a</i>	Study Design Type: Multicenter, historically controlled non-randomized prospective clinical study	Country: United States Follow Up time points: 7-14 days	Study Endpoints: rate of major complications, time to ambulation, time to discharge Performance Outcomes:	Potential Gaps: 6Fr only Diagnostic only Time to hemostasis defined as “the elapsed time between sheath removal and the point hemostasis was first observed” does not align with current definition. Time to ambulation was defined as “the elapsed time between device deployment (or sheath	Justification of Gaps/ How gaps were/are being addressed: The study only included 6Fr per the study design as a pivotal investigation to bring the 6Fr Angio-Seal device to market. The S019 study is a companion study to 8Fr Angio-Seal Early Ambulation in Diagnostic Angiography Patients Study (IDE G940136, PMA 930038/S1) which captures the same data in subjects being treated with the 8Fr undergoing diagnostic

<i>confirmatory study on the STS platform.</i> <i>July 5, 2001 – PMA P930038/S029 FDA approval letter – STS change was approved by FDA – FDA asked for study on information pertinent to the requested change or modification and learned about the device from ongoing or completed studies that may reasonably affect an evaluation of the safety and effectiveness of the device or that may reasonably affect the statement of contraindications, warnings, precautions, and adverse reactions in the draft labeling. Limited to studies sponsored by the applicant</i>				removal in the MP group) and the patients ambulating for 5 minutes or 100 feet” which does not align with current definition. Guidewire included as a device accessory was not evaluated as part of the study. Exclusion criteria: Patients with an access sheath for an invasive procedure that has been in place for more than 36 hours does not align with the current precaution listed in the IFU “ If a patient has had a procedural sheath left in place for longer than 8 hours, consideration should be given to use of prophylactic antibiotics before insertion of the Angio-Seal device. Follow up is 14 days	angiography device compared to manual compression. The data collected regarding the patient population was specific to patients undergoing diagnostic evaluation, however 1 subject was inadvertently included who underwent an interventional procedure. The safety and effectiveness of the Angio-Seal™ VCD has been evaluated since this study in interventional procedures. Additionally, the ongoing PMCF Angio-Seal™ Close Registry enrolls subjects undergoing diagnostic procedures, interventional procedures, or both. To address the definitions of time to hemostasis and time to ambulation (defined as the time, in minutes, between the end of the procedure and the moment when the patient can ambulate, the sponsor has created standardized definitions of time to hemostasis (defined as time, in minutes, between the removal of the procedural sheath up to the first observed arterial/pulsatile bleeding cessation (excluding oozing, which is defined as less than one (1) saturated 10x10cm gauze at access site) in subjects not requiring adjunctive intervention), time to ambulation (ambulation is defined as when the patient is able to stand up by the bedside and walk ten (10) steps unaided, which is equivalent to the patient walking from the hospital bed to the bathroom) and oozing (defined as less than one (1) saturated 10X10cm gauze at access site). These definitions
MDCG 2020-06 Rank: Rank 4 Study Type Per ISO 14155: Confirmatory study	Inclusion Criteria: - At least 18 years of age - Signed informed consent	Patient Demographics: Procedural characteristics:	Safety / Adverse Event Outcomes:		

Study Objective: The primary objective of the study was to determine if the 6F Angio-Seal (A-S) device was safe and effective for the use in closing the arterial puncture following diagnostic angiography (DA). In addition, the study sought to demonstrate the DA patients could be safely ambulated at 30-60 minutes after A-S device deployment and to demonstrate that the rate of complications would not be significantly increased in this group of DA patients who ambulated early compared to a historical manual pressure (MP) treatment control group that ambulated at a later time interval (4-6 hours post sheath removal). Secondary measures also included a comparison of mean time to discharge between the A-S group and a point estimate derived from the historical MP group discharge data.	- Clinically indicated for a diagnostic angiography procedure - Sheath placed in femoral artery for diagnostic procedure is 6F or smaller - Patient should be physically able to walk 100 feet or five minutes immediately following hemostasis procedure. - Patient should be eligible for discharge from the hospital as early as 1-3 hours following device deployment - Patient will be available for follow up office visit 7-14 days after the procedure to be assessed by an investigator or qualified staff. Exclusion criteria - Concomitant femoral venous access in the same leg. - Patient undergoing an interventional cardiac procedure. - Patient is not expected to be physically capable of transferring from the procedure table within ten minutes of sheath removal. - Known allergy to any of the materials used in the device. - Severe acute non-cardiac systemic disease or terminal illness. - Evidence of systemic bacterial or cutaneous infection, including groin infection. - Potential coagulopathy as indicated by platelet count <75,000/pl. - Activated clotting time (ACT) 2 200 seconds (only applicable if patient receives heparin during procedure) Use of thrombolytic agents within 24 hours prior to the	One hundred two (102) patients were enrolled and randomized to the MP control group in the 8F Angio-Seal Early Ambulation Study in 1996 (2 patients were disqualified following randomization and prior to sheath removal). Of the 100 remaining MP patients, 7 had sheaths > 6F; therefore, 93 MP patients were available for comparison to the 6F A-S group. In the 6F Angio-Seal study, 280 patients were enrolled and treated with the 6F A-S device between January 26 and April 29, 1999; no patients were randomized. Of the 280, 80 training patients were not considered part of the pivotal study cohort. Thus, 200 patients were enrolled in the pivotal study, and 6F device placement was attempted. Device placement was abandoned in 2 pivotal study patients due to inability to place the sheath or inability of the device to properly exit the sheath tip. Thus, 198 of these pivotal patients were available for effectiveness analysis (time to hemostasis [TTH], time to ambulation [TIA], or time to discharge [TTD])			will be utilized going forward and are included in the ongoing PMCF Angio-Seal™ Close Registry. This registry also includes data collection on the guidewire that is included with the Angio-Seal™ VIP VCD as a device accessory, which were not collected in the Early Ambulation in Diagnostic Angiography Patients study (S019). The guidewire included as a device accessory was also evaluated in the Clinical Survey Per the PMA submission, the shortened follow window was utilized to capture any complications that may have been missed due to early discharge but would have likely resolved by 14-30 days post procedure. The ongoing PMCF Angio-Seal™ Close Registry includes a 30-day (+ 7 days) post procedure follow up timepoint.

	catheterization procedure that causes the concentration of fibrinogen to be <100mg/dl. - Receiving thrombolytic agents during diagnostic procedure. - Receiving adjunctive platelet inhibitors (i.e. Clopidogrel, Ticlopidine, GP IIb/IIIa inhibitors) or therapeutic anticoagulation therapy (excluding aspirin). - Use of intra-aortic balloon pump support through the puncture site that will be used. Sheath size used for diagnostic procedure is larger than 6F. - Sheath used for diagnostic procedure in place for > 36 hours. - Use of any other investigational drug or device in the last 30 days. - Existence of a palpable hematoma at the end of the diagnostic procedure. - Patient could be pregnant or lactating Post-procedure angiogram demonstrates any of the following contraindications to device placement: • puncture was made in the profunda femoris or at the bifurcation of the common femoral artery. • side-wall entry • posterior wall stick • atherosclerotic vessel • disease in the region of the arteriotomy	compared to the MP data.			

	- vessel tortuosity - femoral artery diameter less than 4 mm - any other condition that could lead to difficulties achieving successful deployment.				

Post-Market Clinical Studies

Many studies have been completed to support the safe and effective use of the Angio-Seal™ VIP.

The table below includes two post-market studies performed with Angio-Seal™ VIP.

Table 2.8: Post Market Confirmatory Studies on Angio-Seal VIP VCD

Post Market Confirmatory Studies on Angio-Seal VIP VCD					
Name of Study / Year:	Study Design Type:	Country:	Study Endpoints: Primary Endpoints:	Potential Gaps	Justification of Gaps/ How gaps were/are being addressed:
Clinical Study of Angio-Seal™ VIP Vascular Closure Device following Diagnostic and/or Interventional Endovascular Procedures (CP0901) (NCT01393041) 2012	Single-center, Post Market Study	Germany <u>Follow Up time points:</u> 30 days post procedure	- Device deployment characteristics & performance - Adverse events (Vascular complications) out to 30-days post-procedure on patients who have undergone a diagnostic and/or interventional procedure - Secondary End-points - Time to hemostasis (TTH) - Rate of minor vascular complications - Performance Outcomes - 100% (49/49) Device deployment successful	6Fr only Time to hemostasis was defined as the time elapsed from device deployment (cutting of the suture) to complete cessation of arterial bleeding, stratified as hemostasis in less than 1 minute, 1-5 minutes, or greater than 5 minutes. No time to ambulation No guidewire <u>Inclusion:</u>	The CP0901 study only captured data on the 6Fr iteration of the device. Since that study, other studies have been conducted on the 8 Fr iteration by itself or have included both the 6Fr and 8Fr iterations of the device. The ongoing PMCF Angio-Seal™ Close Registry collects data on both the 6Fr and 8Fr Angio-Seal™ VIP VCD when used in diagnostic, interventional, or both procedures. The study did not report data on time to ambulation, however Lucatelli

Post Market Confirmatory Studies on Angio-Seal VIP VCD					
			- Hemostasis by device 100% (49/49) - TTH achieved in 1 minute or less in all case	Must be via a retrograde femoral arterial access	2017 and Martin 2008 provide sufficient clinical evidence in support of this clinical benefit. The ongoing PMCF Angio-Seal™ Close Registry collects data on TTH and TTA for the 6Fr and 8Fr Angio-Seal™ VIP VCD when used in diagnostic, interventional, or both procedures. To address the definition of time to hemostasis, the sponsor has created a standardized definitions of time to hemostasis, time to ambulation and oozing. These definitions will be utilized going forward and are included in the ongoing PMCF Angio-Seal™ Close Registry. This registry also includes data collection on time to ambulation and the guidewire that is included with the Angio-Seal™ VIP VCD as a device accessory, which were not collected in the SEAL PM study. The guidewire included as a device accessory was also evaluated in the Clinical Survey “0.035” and 0.038” Guidewire Clinical Survey”
Device Iteration: Angio-Seal™ VIP VCD					Justification of excluding antegrade access. Retrograde is the more common approach. At this juncture in the device development, potential procedural complications may have contaminated the safety endpoint data being captured specifically related to the device.

Post Market Confirmatory Studies on Angio-Seal VIP VCD					
MDCG 2020-06	Inclusion Criteria:	Patient Demographics:	Safety / Adverse Event Outcomes:	Study Objective:	
Rank 2	Patient is undergoing a diagnostic and/or interventional endovascular procedure via a retrograde femoral arterial access	N = 49 patients Mean age 68.5 years Age at baseline (years) 68.5 ± 11.7 (49) (42.1, 89.4) Male 57.1% (28/49) BMI 28.6 kg/m ² History of - Diabetes 42.9% (21/49) - Hypertension 81.6% (40/49) - Acute myocardial infarction 6.1% (3/49) Hypercholesterolemia 34.7% (17/49) - Smoking 28.6% (14/49) - Known coronary artery disease 44.9% (22/49)	- 3 Major vascular complications (6.1%) - 1 Minor vascular complication (1%) that was probably related to the study device	The objective of this study was to assess device deployment performance of the AngioSeal™ and assess adverse events out to 30 days post-procedures on patients undergoing endovascular procedures.	
Study Type Per ISO 14155:	Patient is of legal age		Definitions		
Confirmatory	Patient has given written informed consent for participation prior to the procedure		Major vascular complications: - Vascular injury requiring repair via surgery, angioplasty, - ultrasound-guided compression, - thrombin injection, or other means. - Permanent (unresolved at 30-day post-procedure evaluation) access site-related nerve injury (or access site-related nerve injury requiring intervention); - Access site related bleeding requiring transfusion. - New ipsilateral lower extremity ischemia requiring surgical intervention. - Retroperitoneal bleeding. - Generalized infection requiring prolonged hospitalization and/or treatment with IV antibiotics.	The purpose of this clinical study is to assess device deployment characteristics and performance, as well as adverse events (vascular complications) of the SJM AngioSeal™ VIP 6 French device utilizing the SJM Collagen	
	Patient is willing to undergo all study procedures and adhere to data collection and follow-up requirements				
	Exclusion Criteria:	CVA/TIA 6.1% (3/49) Hypercoagulability conditions 0.0% (0/49) - Bleeding disorders 0.0% (0/49) - Significant peripheral vascular disease including claudication 4.1% (2/49) - Peripheral vascular surgery 0.0% (0/49) - Previous ipsilateral percutaneous femoral access 36.7% (18/49) Arterial puncture closure device used 16.7% (3/18) - Other significant medical history 4.1% (2/49)			
	Patient is unable to provide written informed consent				
	Patient has an inability or is unwilling to adhere to data collection and follow-up requirements				
	Patients who are pregnant or lactating				
	Patient may not participate in another clinical trial which has the potential to impact hemostasis (pharmaceutical/ device/ homeopathic)				
	Patients who initially met the selection criteria and were later determined to be ineligible for the use of the device before, during or following the diagnostic and/ or interventional	Procedural characteristics: Procedure type - Diagnostic 61.2% (30/49) - Interventional 12 (24.5%) - Both 7 (14.3%)	The following are defined as a minor vascular complication: - Access site bleeding requiring ≥ 30 minutes of manual compression to re-achieve hemostasis - Ipsilateral hematoma >10 cm. - Ipsilateral pseudoaneurysm without intervention.		

Post Market Confirmatory Studies on Angio-Seal VIP VCD					
	procedure, due to patient anatomy or circumstances related to the procedure, were designated as "Not Deployed" in the study per the clinical investigational plan. Patients who are designated as "Not Deployed" did not undergo follow-up related to this study. Data from these patients is not pooled with the rest of the clinical data. The primary analysis dataset was based on all patients who signed informed consent, met the selection criteria, and had a VIP device deployed.	Duration of procedure (min) 28.5 ± 22.8 (3, 99) (49) Puncture leg (Right) 93.9% (46/49) Final sheath size: - 6F 100% (49/49) - 6F VIP device 100% (49/49) Preparation/deployment of VIP in accordance with IFU 100% (49/49)	- Ipsilateral arteriovenous fistula. - Ipsilateral deep vein thrombosis. - Local access site infection without prolonged hospitalization.		
Name of Study	Study Design Type:	Country	Study Endpoints:	Potential Gaps:	Justification of Gaps/How gaps were/are being addressed:
Year: SEAL PM Post-Market Study of the St. Jude Medical Angio-Seal™ VIP Vascular Closure Device CV-13-008-EU-VC NCT01858636 2013	Multicentered, Post-Market Study, Characterize Clinical Outcomes of Angio-Seal™ with SJM manufactured collagen	Germany Netherlands Follow Up time points: 30 days	Primary Safety Endpoint: Device or procedure related major vascular complications from implant to 30 days post-procedure Primary Efficacy Endpoint: - Proportion of device deployments that achieve hemostasis within 5 minutes of device deployment Performance Outcomes: - 98.6% (210/213) Devices deployed successfully - 99.1% (210/212) Hemostasis by Device - 0.2 ± 0.3 (210) Time to Hemostasis by device (min) - 98.6% (210/213) Hemostasis within 5 minutes - Device successfully deployed 98.6% (210/213)	Time to hemostasis definition "This is defined as the time elapsed from device deployment (cutting of the suture) to complete cessation of arterial bleeding" exists but does not align with current definition Time to ambulation not defined or collected Guidewire included as a device accessory was not evaluated as part of the study.	There were very limited exclusion criteria as such no potential gaps related to the special patient population (per IFU) identified. To address the definition of time to hemostasis, the sponsor has created a standardized definitions of time to hemostasis, time to ambulation and oozing. These definitions will be utilized going forward and are included in the ongoing PMCF Angio-Seal™ Close Registry. This registry also includes data collection on time to ambulation and the guidewire that is included with the Angio-Seal™ VIP VCD as a device accessory, which were not collected in the SEAL PM study. The study did not report data on time to ambulation, however Lucatelli

Post Market Confirmatory Studies on Angio-Seal VIP VCD					
			- Hemostasis by device 99.1% (210/212) - Time to hemostasis by device (min) 0.2 ± 0.3 (210) (0, 3) - Hemostasis within 5 minutes 98.6% (210/213)		2017 and Martin 2008 provide sufficient clinical evidence in support of this clinical benefit. The ongoing PMCF Angio-Seal™ Close Registry collects data on TTH and TTA for the 6Fr and 8Fr Angio-Seal™ VIP VCD when used in diagnostic, interventional, or both procedures. The guidewire included as a device accessory was also evaluated in the Clinical Survey: 0.035" and 0.038" Guidewire Clinical Survey Report
Device Iteration: Angio-Seal™ VIP VCD					
MDCG 2020-06 Rank: Rank 2	Inclusion Criteria: Subject requires closure of femoral artery puncture resulting from arterial access procedure. Subject is ≥ 18 years of age. Subject is willing to provide written informed consent prior to study device use. Subject is willing and able to adhere to data collection and follow-up requirement. Exclusion Criteria Subject is participating in another clinical trial which has the potential to impact hemostasis. Subject has an active groin infection or systemic infection.	Patient Demographics: 235 Patients Enrolled Age (years) 66.6 ± 11.1 (211) (32, 88) Gender (male) 68.7% (145/211) Height (cm) 172.2 ± 9.1 (211) (146, 205) Weight (kg) 84.2 ± 15.2 (211) (51, 130) Body mass index (kg/m^2) 28.3 ± 4.2 (211) (18.7, 42.2) Body surface area (m^2) 1.97 ± 0.2 (211) (1.44, 2.7) Diabetes 26.1% (55/211) Hypertension 76.3% (161/211) Hypercholesterolemia 59.2% (125/211) Smoking 15.2% (32/211) Known coronary artery disease 55.5% (117/211) Hypercoagulability conditions 26.5% (56/211) Bleeding disorder 0% (0/211)	Safety / Adverse Event Outcomes: - 2.8% (6/211) Patients with Major Vascular Event - Bacterial infections 2/211 (0.9%), 1 device related - VASC bleeding 3/211 (1.4%), 1 device related - VASC pseudoaneurysm 2/211 (0.9%), 2 device-related - Total adjudicated primary endpoints 6/211 (2.8%), 2 device-related Procedural related adverse events - Abdominal pain 1/211 (0.5%) - Angina pectoris 1/211 (0.5%) - Aortic aneurysm 1/211 (0.5%) - Aortic valve stenosis 3/211 (1.4%) - Atrial fibrillation 3/211 (1.4%) - Bacterial infections 3/211 (1.4%)		

Post Market Confirmatory Studies on Angio-Seal VIP VCD					
	Subject has undergone a vascular access procedure on the same leg within the past 90 days. Subject has been previously enrolled in the study. (For subjects undergoing multiple interventions within a timeframe of 90 days, only the initial procedure may constitute a study enrollment.)	Clinically significant peripheral vascular disease, including claudication 4.7% (10/211) Peripheral vascular surgery 0.9% (2/211) Acute myocardial infarction 14.7% (31/211) CVA/TIA 4.7% (10/211) Previous ipsilateral percutaneous femoral access (>90 days prior) 16.6% (35/211) Known complications*2.9% (1/35) Other significant medical history** 55.9% (118/211) <i>* Known complications –the number of procedural complications reported during previous ipsilateral percutaneous femoral access procedures (>90 days prior)</i> <i>** Other frequently reported baseline medical history, as reported by the site, included Atrial Fibrillation, Myocardial Infarction, Renal Insufficiency/Failure, stent procedure, CABG procedure, cardiac arrhythmia, and cardiac valve conditions.</i>	- Cardiomyopathy 1/211 (0.5%) - Coronary artery disease 4/211 (1.9%) - Decreased left ventricular (LV) function 1/211 (0.5%) - Elective cardiac surgery 1/211 (0.5%) - Hematuria 1/211 (0.5%) - Intracerebral hemorrhage 1/211 (0.5%) - Ischemic stroke 1/211 (0.5%) - Myocardial infarction 1/211 (0.5%) - Sick sinus syndrome 1/211 (0.5%) - Sustained ventricular tachycardia 1/211 (0.5%) - Tachycardia 1/211 (0.5%) - Urinary tract infection 1/211 (0.5%) - VASC bleeding 4/211 (1.9%) - VASC hematoma 4/211 (1.9%) - VASC pseudoaneurysm 4/211 (1.9%) - Vascular dissection 1/211 (0.5%)		
		213 Device Deployments	Definitions Major vascular complications Access Site Complications: - Hematoma >10 cm in size requiring surgical or percutaneous intervention - Major bleeding requiring transfusion of ≥2 units of blood or requiring surgical or percutaneous intervention - Pain requiring a hospitalization extended for more than 24 hours or a new hospitalization, or percutaneous or surgical intervention		
		Procedural characteristics: - Diagnostic 72.3% (154/213) Procedure type: Intervention 5.6% (12/213)			

Post Market Confirmatory Studies on Angio-Seal VIP VCD					
		- Procedure type: Diagnostic and Intervention 22.1% (47/213) - Device size: 6F 93% (198/213) - Device size: 8F 7% (15/213) - Puncture leg (Right) 89.2% (190/213) - Puncture leg (Left) 10.8% (23/213) - Procedure duration (min) 28 ± 20.7 (213) (1, 125) - Multiple punctures in same leg 8% (17/213) - Final procedural sheath size: 5F 22.5% (48/213) - Final procedural sheath size: 6F 70.9% (151/213) - Final procedural sheath size: 7F 0.5% (1/213) - Final procedural sheath size: 8F 6.1% (13/213) - ACT prior to sheath removal (sec) 295.3 ± 150.3 (22) (100, 635)	- Infection requiring a hospitalization extended for more than 24 hours or a new hospitalization or treatment with IV antibiotics - A/V Fistula requiring medical intervention (percutaneous or surgical) - Pseudoaneurysm requiring medical intervention (percutaneous or surgical) - Lower Limb Ischemia requiring surgical or medical intervention or resulting in permanent injury/impairments. - Retroperitoneal hemorrhage requiring intervention (percutaneous or surgical)		
Study Type Per ISO 14155: Confirmatory					

Post Market Confirmatory Studies on Angio-Seal VIP VCD		
Study objective The purpose of the study is to characterize the clinical outcomes of Angio-Seal™ with SJM manufactured collagen.		

The table below describes post market animal studies performed on the Angio-Seal™ Evolution device.

Table 2.9: Post-Market Functional Animal Testing

Post-Market Functional Animal Testing					
Name of Study / Year:	Study Design Type:	Study end-points:	Conclusions:	Potential Gaps:	How gaps were/are being addressed:
Animal Study Report, Evolution CER044, Chronic Animal Study 93d, APS #9001467575 2008	Animal study Study objective Document the tissue response and material degradation at early time points for the newest of the Angio-Seal devices and compare to traditional manual closure Species / Additional Details Porcine (2) One access site for each pig was closed using a 6F Angio-Seal™ Evolution device after an exchange of the introducer sheath for a 6F Angio-Seal™ VIP sheath and dilator	Tissue response, Material degradation, Proper vessel patency, optimal healing, vascular inflammation	-There were no adverse changes at any of the sites evaluated: proper vessel patency, optimal healing, and lack of vascular inflammation -The arterial access sites were healed with organized and mature neointima at the puncture site -There were no appreciable differences between the control and test sites at 93 days (no to minimal residual perivascular hemosiderin was present and fibrosis was minimal to mild) -Implantation of the St. Jude Medical (SJM) Angio-Seal Evolution Vascular Closure Device at a vascular access site followed by a 93-day observation period showed no adverse changes and no appreciable qualitative differences between treated and manually compressed vessels	- Study utilized only 6Fr devices - Animal model that is single species - Does not use the entire Angio-Seal VIP device Justification of Gaps - Preclinical testing of newly developed vascular in situ tissue material testing in animal models is commonplace - Anticipated clinical difference between 6Fr and 8Fr are negligible.	Animal Study at APS #90017642 &, Acute Animal Study at APS #90018838 include 8 Fr devices.

<u>Device Iteration:</u> Angio-Seal™ Evolution Angio-Seal™ VIP sheath and dilator <u>MDCG 2020-06 Rank:</u> Rank 11, Animal Study <u>Study Type Per ISO 14155:</u> <u>N/A</u>					
<u>Name of Study / Year:</u> Animal Study Report, Evolution CER056, Acute Animal Study at APS #90017642 2008	<u>Study Design Type:</u> Animal study <u>Study objective</u> Evaluate collagen placement relative to the control groups as well as deployment observations, such as forces to deploy the devices <u>Species / Additional Details</u> Porcine (6)	<u>Study endpoints:</u> Evaluation of collagen placement, device deployment observations (forces to deploy device), Hemostasis, Time to hemostasis	<u>Conclusions:</u> - No suture release mechanism activation occurred on test or control devices in this study. One FW control device was noted to have high force, but the suture release mechanism did not activate - All 6 FW devices had collagen that resulted in final placement above the skin and had a dimpling effect. All NFW and VIP devices had no collagen visible above skin level - 2 of 12 NFW devices required 4 and 5 minutes of manual compression. All FW devices had immediate hemostasis. One VIP device required 2.5 minutes of manual compression - Average collagen height for the FW devices placed in the artery was 10.03mm. VIP average collagen height was 7.17mm and NFW devices averaged 6.96mm. This is the	<u>Potential Gaps:</u> - Definition of ‘hemostasis’ in the study is ‘the cessation of blood flow from the arteriotomy.’ - Definition of ‘time to hemostasis’ is the time from the device deployment to the time hemostasis is achieved.” - Single species animal model – Porcine arteries are approximately 0.5mm thick, half of the thickness of human arteries. <u>Justification of Gaps</u> 1. The porcine arteries are very challenging since they are approximately 0.5mm thick which is about half the thickness of human arteries. The porcine model is an aggressive model for closure devices because the tissue is more muscular than a human and the arteries will tear easier than human arteries. The femoral arteries also dive down so the angle of	<u>How gaps were/are being addressed:</u> 1. Not clinical data per XXX and do not need to address this gap

			height of the collagen from the top of the artery to the knot - During cut down it was noted that all collagen was sandwiched down directly on top of the arterial wall, with the exception of one 6F NFW that had muscular tissue (measured at 2.8mm) between the collagen and the artery resulting in a Hematoma, and one 6F FW device had tissue between the anchor and the collagen, but a measurement was not recorded - There were 3 hematomas noted during cut down to the arteries. These were noted in 1 6F and 2 8F NFW device groups - 4 artery tears were noted in the 8F NFW group, two of which were Hi Ho deployments (Anchor half in and half out of the arterial wall)	sheath insertion is often shallow. In a Hi Ho deployment half of the anchor is inside the vessel and half is above due to the tear.	
<u>Device Iteration:</u> Angio-Seal™ Evolution / P1076 collagen; Angio-Seal™ VIP <u>MDCG 2020-06 Rank:</u> Rank 11, Animal Study <u>Study Type Per ISO 14155:</u> <u>N/A</u>					
<u>Name of Study / Year:</u> Animal Study Report, Evolution CER056, Acute Animal	<u>Study Design Type:</u> Animal study <u>Study objective</u>	<u>Study end-points:</u> Evaluation of collagen placement, device deployment	<u>Conclusions:</u> - 8 total deployments were completed in 2 animals.	<u>Potential Gaps:</u> 1. Only on 8Fr devices 2. Definition of 'hemostasis' in the study is 'the cessation of blood	<u>How gaps were/are being addressed:</u> Animal Study Report,

Study at APS #90018838 2008	Evaluate the Evolution Vascular Closure device in the porcine model and to confirm previously acquired bench model results in the in-vitro tissue model. We compared the Evolution device using P1076 collagen and the Angio Seal VIP device Hemostasis and Time to hemostasis <u>Species / Additional Details</u> Porcine (2/ Up to 4 deployments per animal) 4 Evolution 8F test devices were placed along with 4 VIP 8F controls	observations (forces to deploy device), Hemostasis, Time to hemostasis	- 4 Evolution 8F test devices were placed along with 4 VIP 8F controls - 2 of 4 VIP devices required manual compression - During cut down it was noted that all collagen was sandwiched down directly on top of the arterial wall, with the exception of one 8F VIP that had muscular tissue (measured at 2-3mm) between the collagen and the artery resulting in the need for manual compression - There was 1 hematoma noted during cut down to the arteries. This occurred in the 8F Evolution device which was deployed subcutaneously in the tissue - No Hi Ho deployments were observed in the test or control group. 2 artery tears were noted in the 8F EVO group. 3 artery tears were noted in the 8F VIP group	flow from the arteriotomy.’ 3. Definition of ‘time to hemostasis’ is the time from the device deployment to the time hemostasis is achieved.” 4. Single species animal model – Porcine arteries are approximately 0.5mm thick, half of the thickness of human arteries. <u>Justification of Gaps</u> -The porcine arteries are very challenging since they are approximately 0.5mm thick which is about half the thickness of human arteries. The porcine model is an aggressive model for closure devices because the tissue is more muscular than a human and the arteries will tear easier than human arteries. The femoral arteries also dive down so the angle of sheath insertion is often shallow. In a Hi Ho deployment half of the anchor is inside the vessel and half is above due to the tear.	Evolution CER044, Chronic Animal Study 93d, APS #9001467575 & Animal Study Report, Evolution CER056, Acute Animal Study at APS #90017642 include 6Fr devices
<u>Device Iteration:</u> Angio-Seal™ Evolution / P1076 collagen; Angio-Seal™ VIP <u>MDCG 2020-06 Rank:</u> Rank 11, Animal Study <u>Study Type Per ISO 14155:</u> <u>N/A</u>					

Post-Market Clinical Follow Up

Terumo Medical Corporation continuously collects information on the device after it is released in market to demonstrate long-term safety and performance of the device. Terumo is currently conducting a post-market clinical investigation. There were two high quality clinical surveys performed to document feedback related to the use of Angio-Seal™ VIP.

The table below summarizes Post Market Clinical Follow Up Activities on Angio-Seal™ VIP VCD:

Table 2.10: Post Market Clinical Follow Up Activities on Angio-Seal™ VIP VCD

Post Market Clinical Follow Up Activities on Angio-Seal VIP VCD			
Name of Activity	Study Design Type:	Type of procedures	1-5 Scale device evaluations – Survey statements
20210293 – AngioSeal™ VIP Vascular Closure System Clinical Survey Report (2020)	<u>Study objective</u> The purpose of this report is to document prospective feedback related to the use of Angio-Seal™ VIP Vascular Closure Device (AS VIP VCD) distributed in the European Union's market. This was done in the form of a clinical survey and completed by physicians after they used the device in a case procedure. The feedback includes technical outcomes, clinical performance, and safety pertaining to the performance of the device.	- Coronary - Peripheral - Visceral - Peripheral/coronary - Visceral/Peripheral - Coronary/other - Other (Neuro, TAVR, Stroke, Stroke thrombectomy, carotid stent, AVM, neurointerventional, profunda branch bleed, structural, and	<u>Acceptance criteria</u> • Ratings 3-5 are deemed favorable • All statements should have ≥ 80% positive feedback
<u>Device Iteration:</u> Angio-Seal™ VIP VCD			<u>Outcomes</u> • No survey statement received > 20% unfavorable responses No complaints against product were filled
<u>MDCG 2020-06 Rank:</u> Rank 4 – High quality survey			<u>Results</u> Although 138 surveys were completed in total, as discussed in Deviation 4 of Section 4, six surveys were excluded and thus, a total of 132 surveys (90 6Fr Angio-Seal™ VIP VCD and 42 8Fr Angio-Seal™ VIP VCD) were considered valid for this evaluation.
<u>Study Type Per ISO 14155:</u> N/A	This device evaluation is performed to support AS VIP VCD European Medical Device Regulation (EU MDR 2017/745) submission to prove continued safe and effective intended use of the product		<u>Conclusions</u> No survey statement received more than 20% unfavorable responses and are deemed acceptable per the protocol. Thus, based on the acceptance criteria defined in the protocol and physician feedback based on the ratings and comments recorded during the surveys, the Angio-Seal™ VIP Vascular Closure Device has successfully demonstrated that it

Post Market Clinical Follow Up Activities on Angio-Seal VIP VCD			
	and further document sufficient clinical evidence for AS VIP VCD per MDCG 2020-6 (Guidance on sufficient clinical evidence for legacy devices). High Quality clinical surveys are identified in appendix III of MDCG 2020-6 as an additional method of demonstrating sufficient clinical evidence for a legacy device under EU MDR 2017/745. 138 cases 88 unique physicians with > 5 years experiences 95 surveys completed for 6 Fr Angio-Seal™ VIP, and 43 completed for 8 Fr Angio-Seal™ VIP VCD	cerebral diagnostic) Devices included - 6Fr Angio-Seal™ VIP VCD OUS / 610132 - 8Fr Angio-Seal™ VIP VCD OUS / 610133 - 6Fr Angio-Seal™ VIP VCD US / 610130 - 8Fr Angio-Seal™ VIP VCD US / 610131 Participating Physician Characteristics - Endovascular and operative neurovascular surgeon - Interventional cardiologist - Interventional neurologists - Interventional neuroradiologists - Interventional radiologist - Neurosurgeon - Vascular surgeon - Neurointerventional radiologist	continues to allow for safe and effective intended use

Table 2.11 Planned/Ongoing Post Market Clinical Follow Up Activities on Angio-Seal VIP

Name of Study /	Study Design Type:	Country	Study endpoints:	Potential Gaps	How gaps were/are being addressed:
Year: Angio-Seal™ Close Post-Market Clinical Investigation of the Angio-Seal™ VIP VCD: A Prospective, Multi-Center, Observational Study in Europe NCT0533525 CIP # T139E4 Ongoing/ Enrolling subjects Device Iteration: Angio-Seal™ VIP VCD	Prospective, Multi-center, Observational, Post-market observational study	Germany, France, Netherlands Follow Up time points: 30 + 7 days post procedure Duration First patient enrolled: September 2022 Last patient enrolled: December 2022 Last patient last contact: July 2023	Primary endpoints: Composite safety and effectiveness endpoint of successful Angio-Seal™ VIP VCD puncture site hemostasis and freedom from major complications of the access site limb within 6 hours post-procedure, where: Successful puncture site hemostasis is defined as: -Cessation of arterial bleeding (excluding oozing) achieved in the interventional lab in subjects not requiring adjunctive intervention at the access site, including manual compression. Major complications attributable to Angio-Seal™ VIP VCD are defined as: -Access site-related bleeding (BARC type 2, 3 or 5, following Angio-Seal™ VIP VCD deployment, as per Bleeding Academic Research Consortium (BARC) classification) -Femoral puncture site hematoma >6 -Pseudoaneurysms requiring intervention -Femoral access site arteriovenous fistulas -Access site infection requiring hospitalization -Embolism (due to Anchor fracture) -Thrombosis at puncture site (due to collagen plug disposition into the artery) -Allergic Reaction to Angio-Seal™ VIP VCD components -Foreign body reaction -Inflammation and Edema Secondary endpoints 1. Freedom from any minor complications at the target limb access site within 6 hours post-procedure. Minor complications are defined as: -Femoral puncture site Hematoma <6cm -Access site Infection not requiring hospitalization		

			-Pseudoaneurysm not requiring intervention - Vasovagal response 2. Freedom from any major and minor complications at the target limb access site after 6 hours to 30 days post-procedure. 3. Time to hemostasis (TTH): defined as time, in minutes, between the removal of the procedural sheath up to the first observed arterial/pulsatile bleeding cessation (excluding oozing) in subjects not requiring adjunctive intervention in the interventional lab 4. Time to ambulation (TTA): defined as the time, in minutes, between the end of the procedure and the moment when the patient is able to ambulate (ambulation is defined as when the patient is able to stand up by the bedside and walk ten (10) steps unaided, which is equivalent to the patient walking from the hospital bed to the bathroom) 5. Quality of Life assessment (EQ-5D) 6. Angio-Seal™ VIP VCD usability <u>Performance Outcomes</u>		
MDCG 2020-06 Rank: Rank 3 Study Type Per ISO 14155: Exploratory (Angio-Seal™ VIP VCD)	Inclusion Criteria: 1) Subject is ≥ 18 years old 2) Subject is willing and able to complete the follow-up requirements 3) Subject has the mental capacity (i.e., does not require the use of a Legally Authorized Representative) to sign the study Informed Consent Form (ICF)	Patient Demographics: <i>To be described upon completion of the registry</i> Procedural Characteristics: <i>To be described upon completion of the registry</i>	Safety / Adverse Event Outcomes:	Justification of Gaps:	

Study objective: To evaluate the safety and effectiveness of the Angio-Seal™ VIP VCD in patients undergoing endovascular procedures via femoral access in real-world setting. Angio-Seal™ VIP VCD is a medical device indicated for use in closing and reducing time to hemostasis at the femoral arterial puncture site in patients who have undergone diagnostic angiography procedures or interventional procedures (minimal artery diameter of 4mm).	4) Patient undergoing diagnostic or interventional endovascular procedure compatible with the use of Angio-Seal™ VIP VCD 5) Puncture site located at the femoral artery (i.e., between the inguinal ligament and the bifurcation of the superficial femoral and profunda femoris arteries) 6) Angio-Seal™ VIP VCD deployed as per instruction for use by a trained operator Exclusion Criteria: Use of the Angio-Seal™ VIP VCD on puncture sites other than the femoral artery • Repuncture of the femoral artery within 90 days at the same access site • Lumen diameter of femoral artery < 4 mm • Patients with clinically significant peripheral vascular disease at the puncture site (luminal narrowing of >40% within 5 mm of the puncture site) • Puncture site at or distal to the bifurcation of the superficial femoral and profunda femoris artery • Puncture site proximal to the inguinal ligament • Procedure sheath placed through the superficial femoral artery into the profunda femoris • Multiple femoral punctures • Known or suspected posterior femoral wall puncture • Use of > 8F primary introducer sheaths or devices for an 8Fr Angio-Seal™ VIP VCD, or use of a >6Fr primary introducer sheaths or devices for a 6Fr Angio-Seal™ VIP VCD • Any condition that would make use of Angio-Seal™ VIP VCD inappropriate	Planning on 230 enrolled subjects undergoing diagnostic or interventional endovascular procedures in which Angio-Seal™ VIP VCD is used			

	(as per IFU and investigators' discretion)				

1.4.3 Clinical Background of the Device

The Angio-Seal™ VIP VCD has a proven clinical track record of safety and performance. The Angio-Seal™ VIP VCD has been on the market for over 17 years and is considered lacking in novelty.

1.4.4 Clinical Evidence Supporting CE Marking

Angio-Seal™ approval is based on pre-clinical and clinical studies. Studies have shown:

- Angio-Seal™ VIP VCD reduces the time to stop the bleeding (hemostasis).
- Angio-Seal™ VIP VCD reduces the time to regain the ability to walk (ambulation). Unwanted side-effects happen about 5.2% of the time.
- 5,000 patients experienced unwanted side-effects.

1.4.5 Safety and Performance

Angio-Seal™ VIP VCD has a shorter time to stop the bleeding (hemostasis) at the procedure site compared to manual pressure. Angio-Seal™ VIP VCD reduces the time to regain the ability to walk (ambulation) compared to manual pressure. Angio-Seal™ VIP VCD and manual pressure have similar risks. The clinical benefits of using Angio-Seal™ VIP VCD outweigh the risks. Your doctor has considered the risks and decided that this device is appropriate for your procedure. See Section 2.1 for patient population information. If you have an existing health condition, speak with your doctor.

There were 2 post-market clinical studies for Angio-Seal™. The 2 studies included a 2012 study with 49 patients and a 2013 study with 235 patients. These studies confirmed the safety and effectiveness of Angio-Seal™. These studies also confirmed that Angio-Seal™ performs as it is intended.

Terumo Medical Corporation continues to collect post-market data on the device. This is to demonstrate long-term safety and performance of the device. Terumo is performing a post-market clinical follow up study. This study will gather additional information for safety and performance of the Angio-Seal™ device. In addition, a survey of physicians will be performed to gather information on the guidewire used with the device unit. The survey will gather information from physicians on what type of guidewire they use during the procedure.

1.5 Possible Therapeutic Alternatives

Manual pressure is the standard method to stop bleeding. It is the use of external hand pressure, by a member of the hospital staff, to the hole in the artery. There are also mechanical devices to stop bleeding. An alternative therapy, such as manual or mechanical compression, would be used when the conditions for safe deployment of Angio-Seal™ VIP VCD are not met.

When considering other treatments, contact your doctor who can consider your medical history.

1.6 Suggested User Profile and Training

Your doctor has been trained by Terumo Medical Corporation for the use of this device (Angio-Seal™ program).

2.0 BIBLIOGRAPHY

Clinical Background Publications

1. Arunothayaraj S, Palmer S, Whitbourn R, et al. The Art of the Seal. J Am Coll Cardiol Interv. 2021 Jun, 14 (11) e117–e118. <https://doi.org/10.1016/j.jcin.2020.12.043>
2. E Kim, BG Sebastiao, A Lee, S Ande, J Shankar. Safety and effectiveness of vascular closure devices in interventional radiological procedures. Interventional neuroradiology: journal of peritherapeutic neuroradiology, surgical procedures and related neurosciences. 2022. :#pages#
3. F Frenzel, P Fries, R Shayesteh-Kheslat, A Buecker, A Massmann. Single Angio-Seal vascular closure device for transfemoral access exceeding 8F. Journal of cardiology. 2020. 76:#pages#
4. G Zenunaj, L Traina, P Acciarri, M Mucignat, S Scian, F Alesiani, R Serra, V Gasbarro. Superficial Femoral Artery Access for Infrainguinal Antegrade Endovascular Interventions in the Hostile Groin: A Prospective Randomized Study. Annals of vascular surgery. 2022. :#pages#
5. Kreutz, Rolf P., et al. "Femoral artery closure devices vs manual compression during cardiac catheterization and percutaneous coronary intervention." Journal of the Society for Cardiovascular Angiography & Interventions (2022): 100370.
6. LE Jones, KH Yang, RW Feldtman, PV Uceda, CA Ferrara, JM Caruso, JL Richmond, SS Ahn. Safety and Efficacy of Arterial Closure Devices in an Office-Based Angiosuite. Annals of vascular surgery. 2018. 51:#pages#
7. N Ito, K Hishikari, H Hikita, A Takahashi, T Yonetsu, T Sasano. Images of angioscopy and intravascular ultrasound for dislodged Angio-Seal VIP vascular closure device. Cardiovascular intervention and therapeutics. 2022. 37:#pages#
8. M. Shukor, et al. "Managing Unintentional Iatrogenic Triple Lumen Central Venous Catheter (TLCVC) Arterial Puncture."The Interventionalist Journal 1.2 (2021): 12-18.
9. PS Dhillon, N McConachie, R Lenthall, S Nair, W Izzath. Use of 8Fr angio-seal for closure of femoral arteriotomy following use of 8Fr and 9Fr sheaths in patients undergoing mechanical thrombectomy for acute ischaemic stroke. Journal of clinical neuroscience : official journal of the Neurosurgical Society of Australasia. 2021. 88:#pages#
10. Shah, Aakash, et al. "Complete percutaneous decannulation from femoral venoarterial extracorporeal membrane oxygenation. "JTCVS techniques 6 (2021): 75-81.

11. Woźniak, Witold, and Piotr Ciostek. "Infectious arteritis of the femoral artery as a cause of recurrent hemorrhage after endovascular treatment of thrombosed popliteal artery aneurysm." *Acta Angiologica* 26.3 (2020): 113-117.
12. Su SF, Chang MY, Wu MS, Liao YC. Safety and efficacy of using vascular closure devices for hemostasis on sheath removal after a transfemoral artery percutaneous coronary intervention. *Jpn J Nurs Sci*. 2019 Apr;16(2):172-183. doi: 10.1111/jjns.12221. Epub 2018 Jul 25.

3.0 REVISION HISTORY

SSCP revision number	Date issued	Change description	Revision validated by the Notified Body
SSCP-AS-2022 – Revision 1	May 9, 2023	Initial Issue, upon MDR CE certification	☒ Yes Validation language: English ☐ No (only applicable for class IIa or some IIb implantable devices (MDR, Article 52 (4) 2nd paragraph) for which the SSCP is not yet validated by the NB) ☐ Other (explain):
SSCP-AS-2023 – Revision 1	February 15, 2024	Updated for the period between January 1, 2018 through April 30, 2023	☐ Yes Validation language: English ☐ No (only applicable for class IIa or some IIb implantable devices (MDR, Article 52 (4) 2nd paragraph) for which the SSCP is not yet validated by the NB) ☒ Other (explain): This version was updated and not required to be submitted to NSAI per their communication. The next version (SSCP-AS-2024) will be submitted to NSAI and validated.

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1. Arunothayaraj S, Palmer S, Whitbourn R, et al. The Art of the Seal. J Am Coll Cardiol Interv. 2021 Jun, 14 (11) e117–e118. <https://doi.org/10.1016/j.jcin.2020.12.043>
2. E Kim, BG Sebastiao, A Lee, S Ande, J Shankar. Safety and effectiveness of vascular closure devices in interventional radiological procedures. Interventional neuroradiology: journal of peritherapeutic neuroradiology, surgical procedures and related neurosciences. 2022. :#pages#
3. F Frenzel, P Fries, R Shayesteh-Kheslat, A Buecker, A Massmann. Single Angio-Seal vascular closure device for transfemoral access exceeding 8F. Journal of cardiology. 2020. 76:#pages#
4. G Zenunaj, L Traina, P Acciarri, M Mucignat, S Scian, F Alesiani, R Serra, V Gasbarro. Superficial Femoral Artery Access for Infrainguinal Antegrade Endovascular Interventions in the Hostile Groin: A Prospective Randomized Study. Annals of vascular surgery. 2022. :#pages#
5. Kreutz, Rolf P., et al. "Femoral artery closure devices vs manual compression during cardiac catheterization and percutaneous coronary intervention." Journal of the Society for Cardiovascular Angiography & Interventions (2022): 100370.
6. LE Jones, KH Yang, RW Feldtman, PV Uceda, CA Ferrara, JM Caruso, JL Richmond, SS Ahn. Safety and Efficacy of Arterial Closure Devices in an Office-Based Angiosuite. Annals of vascular surgery. 2018. 51:#pages#
7. N Ito, K Hishikari, H Hikita, A Takahashi, T Yonetsu, T Sasano. Images of angioscopy and intravascular ultrasound for dislodged Angio-Seal VIP vascular closure device. Cardiovascular intervention and therapeutics. 2022. 37:#pages#
8. M. Shukor, et al. "Managing Unintentional Iatrogenic Triple Lumen Central Venous Catheter (TLCVC) Arterial Puncture." The Interventionalist Journal 1.2 (2021): 12-18.
9. PS Dhillon, N McConachie, R Lenthall, S Nair, W Izzath. Use of 8Fr angio-seal for closure of femoral arteriotomy following use of 8Fr and 9Fr sheaths in patients undergoing mechanical thrombectomy for acute ischaemic stroke. Journal of clinical neuroscience : official journal of the Neurosurgical Society of Australasia. 2021. 88:#pages#
10. Shah, Aakash, et al. "Complete percutaneous decannulation from femoral venoarterial extracorporeal membrane oxygenation." JTCVS techniques 6 (2021): 75-81.

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2.0 BIBLIOGRAFI

Kliniske baggrundspublikationer

1. Arunothayaraj S, Palmer S, Whitbourn R, et al. The Art of the Seal. J Am Coll Cardiol Interv. 2021 Jun, 14 (11) e117–e118. <https://doi.org/10.1016/j.jcin.2020.12.043>
2. E Kim, BG Sebastiao, A Lee, S Ande, J Shankar. Safety and effectiveness of vascular closure devices in interventional radiological procedures. Interventional neuroradiology: journal of peritherapeutic neuroradiology, surgical procedures and related neurosciences. 2022. :#pages#
3. F Frenzel, P Fries, R Shayesteh-Kheslat, A Buecker, A Massmann. Single Angio-Seal vascular closure device for transfemoral access exceeding 8F. Journal of cardiology. 2020. 76:#pages#
4. G Zenunaj, L Traina, P Acciarri, M Mucignat, S Scian, F Alesiani, R Serra, V Gasbarro. Superficial Femoral Artery Access for Infringuinal Antegrade Endovascular Interventions in the Hostile Groin: A Prospective Randomized Study. Annals of vascular surgery. 2022. :#pages#
5. Kreutz, Rolf P., et al. "Femoral artery closure devices vs manual compression during cardiac catheterization and percutaneous coronary intervention." Journal of the Society for Cardiovascular Angiography & Interventions (2022): 100370.
6. LE Jones, KH Yang, RW Feldtman, PV Uceda, CA Ferrara, JM Caruso, JL Richmond, SS Ahn. Safety and Efficacy of Arterial Closure Devices in an Office-Based Angiosuite. Annals of vascular surgery. 2018. 51:#pages#
7. N Ito, K Hishikari, H Hikita, A Takahashi, T Yonetsu, T Sasano. Images of angioscopy and intravascular ultrasound for dislodged Angio-Seal VIP vascular closure device. Cardiovascular intervention and therapeutics. 2022. 37:#pages#
8. M. Shukor, et al. "Managing Unintentional Iatrogenic Triple Lumen Central Venous Catheter (TLCVC) Arterial Puncture."The Interventionalist Journal 1.2 (2021): 12-18.
9. PS Dhillon, N McConachie, R Lenthall, S Nair, W Izzath. Use of 8Fr angio-seal for closure of femoral arteriotomy following use of 8Fr and 9Fr sheaths in patients undergoing mechanical thrombectomy for acute ischaemic stroke. Journal of clinical neuroscience : official journal of the Neurosurgical Society of Australasia. 2021. 88:#pages#
10. Shah, Aakash, et al. "Complete percutaneous decannulation from femoral venoarterial extracorporeal membrane oxygenation. "JTCVS techniques 6 (2021): 75-81.

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Première étude clinique exploratoire du Angio-Seal (HPCD) - itération initiale pré-commercialisation chez des humains (6 French)	Itération du Angio-Seal	MDCG 2020-06 Rank
6 French Pilot Study	Angio-Seal Hemostatic Closure Devices (itération initiale)	Rank 4
Études de confirmation précommercialisation – itération initiale (HPCD)	Itération du Angio-Seal	MDCG 2020-06 Rank
US Phase II Europe Phase II Early Ambulation in Diagnostic Angiography Patients	Angio-Seal Hemostatic Closure Devices (itération initiale)	Rank 4
Études d’approbation précommercialisation (menant au Angio-Seal VIP VCD)	Itération du Angio-Seal	MDCG 2020-06 Rank
Clinical Evaluation of the Angio-Seal Self-Tightening Slipknot (STS) Hemostatic Puncture Closure Device (IDE G990306)	Angio-Seal STS	Rank 2
Angio-Seal STS Early Ambulation and Discharge Study	Angio-Seal STS	Rank 2
Études cliniques post-commercialisation terminées	Itération du Angio-Seal	MDCG 2020-06 Rank
Clinical Study of Angio-Seal™ VIP Vascular Closure Device following Diagnostic and/or Interventional Endovascular Procedures (CP0901) (2009)	Angio-Seal™ VIP VCD	Rank 2
SEAL PM Post-Market Study of the St. Jude Medical Angio-Seal™ VIP Vascular Closure Device CV-13-008-EU-VC (2013)	Angio-Seal™ VIP VCD	Rank 2
Activités de suivi clinique post-commercialisation terminées	Itération du Angio-Seal	MDCG 2020-06 Rank
20210293 – Angio-Seal VIP Vascular Closure Device Clinical Survey Report (2020)	Angio-Seal™ VIP VCD	Rank 4
20210298 – Guidewire Clinical Survey Report (2020)	Angio-Seal™ VIP VCD	Rank 4
Essais post-commercialisation sur des animaux	Itération du Angio-Seal	MDCG 2020-06 Rank
Angio-Seal® Vascular Closure Device Biological Evaluation & Testing Report 17-3-121 (2018)	Angio-Seal™ VIP VCD Angio-Seal™ Evolution	Rank 12
Animal Study Report, Evolution CER044, Chronic Animal Study 93d, APS #90014675 (2008)	Angio-Seal™ VIP VCD Angio-Seal™ Evolution	Rank 11
Animal Study Report, Evolution CER056, Acute Animal Study at APS #90017642 (2008)	Angio-Seal™ VIP VCD Angio-Seal™ Evolution	Rank 11
Animal Study Report, Evolution CER056, Acute Animal Study at APS #90018838 (2008)	Angio-Seal™ VIP VCD Angio-Seal™ Evolution	Rank 11
Études cliniques post-commercialisation planifiées/en cours	Itération du Angio-Seal	MDCG 2020-06 Rank
The Angio-Seal Close: Post-market Clinical Investigation of the Angio-Seal™ VIP Vascular Closure Device (VCD): A Prospective, Multi-Center, Observational Study in Europe (T139E4)	Angio-Seal™ VIP VCD	Rank 3 (ongoing)

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Završena klinička ispitivanja nakon stavljanja u promet	Iteracija proizvoda Angio-Seal	MDCG 2020-06 Rank
Clinical Study of Angio-Seal™ VIP Vascular Closure Device following Diagnostic and/or Interventional Endovascular Procedures (CP0901) (2009)	Angio-Seal™ VIP VCD	Rank 2
SEAL PM Post-Market Study of the St. Jude Medical Angio-Seal™ VIP Vascular Closure Device CV-13-008-EU-VC (2013)	Angio-Seal™ VIP VCD	Rank 2
Obavljene aktivnosti kliničkog praćenja nakon stavljanja u promet	Iteracija proizvoda Angio-Seal	MDCG 2020-06 Rank
20210293 – Angio-Seal VIP Vascular Closure Device Clinical Survey Report (2020)	Angio-Seal™ VIP VCD	Rank 4
20210298 – Guidewire Clinical Survey Report (2020)	Angio-Seal™ VIP VCD	Rank 4
Testiranje na životinjama nakon stavljanja u promet	Iteracija proizvoda Angio-Seal	MDCG 2020-06 Rank
Angio-Seal® Vascular Closure Device Biological Evaluation & Testing Report 17-3-121 (2018)	Angio-Seal™ VIP VCD Angio-Seal™ Evolution	Rank 12
Animal Study Report, Evolution CER044, Chronic Animal Study 93d, APS #90014675 (2008)	Angio-Seal™ VIP VCD Angio-Seal™ Evolution	Rank 11
Animal Study Report, Evolution CER056, Acute Animal Study at APS #90017642 (2008)	Angio-Seal™ VIP VCD Angio-Seal™ Evolution	Rank 11
Animal Study Report, Evolution CER056, Acute Animal Study at APS #90018838 (2008)	Angio-Seal™ VIP VCD Angio-Seal™ Evolution	Rank 11
Planirana/tekuća klinička ispitivanja nakon stavljanja u promet	Iteracija proizvoda Angio-Seal	MDCG 2020-06 Rank
The Angio-Seal Close: Post-market Clinical Investigation of the Angio-Seal™ VIP Vascular Closure Device (VCD): A Prospective, Multi-Center, Observational Study in Europe (T139E4)	Angio-Seal™ VIP VCD	Rank 3 (ongoing)

Angio-Seal™ VIP pripada liniji medicinskih proizvoda Angio-Seal™ koji su od 1994. godine prisutni u Europi te od 1996. u SAD-u. Cilj kliničkih ispitivanja bio je potvrditi sigurnost i učinkovitost proizvoda Angio-Seal™.

Tablica 2.5: Prva ispitivanja na ljudima i ispitivanja prije stavljanja u promet

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2.0 IRODALOMJEGYZÉK

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1. Arunothayaraj S, Palmer S, Whitbourn R, et al. The Art of the Seal. J Am Coll Cardiol Intv. 2021 Jun, 14 (11) e117–e118. <https://doi.org/10.1016/j.jcin.2020.12.043>
2. E Kim, BG Sebastiao, A Lee, S Ande, J Shankar. Safety and effectiveness of vascular closure devices in interventional radiological procedures. Interventional neuroradiology: journal of peritherapeutic neuroradiology, surgical procedures and related neurosciences. 2022. :#pages#
3. F Frenzel, P Fries, R Shayesteh-Kheslat, A Buecker, A Massmann. Single Angio-Seal vascular closure device for transfemoral access exceeding 8F. Journal of cardiology. 2020. 76:#pages#
4. G Zenunaj, L Traina, P Acciarri, M Mucignat, S Scian, F Alesiani, R Serra, V Gasbarro. Superficial Femoral Artery Access for Infrainguinal Antegrade Endovascular Interventions in the Hostile Groin: A Prospective Randomized Study. Annals of vascular surgery. 2022. :#pages#
5. Kreutz, Rolf P., et al. "Femoral artery closure devices vs manual compression during cardiac catheterization and percutaneous coronary intervention." Journal of the Society for Cardiovascular Angiography & Interventions (2022): 100370.
6. LE Jones, KH Yang, RW Feldtman, PV Uceda, CA Ferrara, JM Caruso, JL Richmond, SS Ahn. Safety and Efficacy of Arterial Closure Devices in an Office-Based Angiosuite. Annals of vascular surgery. 2018. 51:#pages#
7. N Ito, K Hishikari, H Hikita, A Takahashi, T Yonetsu, T Sasano. Images of angioscopy and intravascular ultrasound for dislodged Angio-Seal VIP vascular closure device. Cardiovascular intervention and therapeutics. 2022. 37:#pages#
8. M. Shukor, et al. "Managing Unintentional Iatrogenic Triple Lumen Central Venous Catheter (TLCVC) Arterial Puncture."The Interventionalist Journal 1.2 (2021): 12-18.
9. PS Dhillon, N McConachie, R Lenthall, S Nair, W Izzath. Use of 8Fr angio-seal for closure of femoral arteriotomy following use of 8Fr and 9Fr sheaths in patients undergoing mechanical thrombectomy for acute ischaemic stroke. Journal of clinical neuroscience : official journal of the Neurosurgical Society of Australasia. 2021. 88:#pages#
10. Shah, Aakash, et al. "Complete percutaneous decannulation from femoral venoarterial extracorporeal membrane oxygenation. "JTCVS techniques 6 (2021): 75-81.

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Studi clinici post-commercializzazione completati	Iteration of Angio-Seal	MDCG 2020-06 Rank
Clinical Study of Angio-Seal™ VIP Vascular Closure Device following Diagnostic and/or Interventional Endovascular Procedures (CP0901) (2009)	Angio-Seal™ VIP VCD	Rank 2
SEAL PM Post-Market Study of the St. Jude Medical Angio-Seal™ VIP Vascular Closure Device CV-13-008-EU-VC (2013)	Angio-Seal™ VIP VCD	Rank 2
Attività di follow-up clinico post-commercializzazione completate	Iteration of Angio-Seal	MDCG 2020-06 Rank
20210293 – Angio-Seal VIP Vascular Closure Device Clinical Survey Report (2020)	Angio-Seal™ VIP VCD	Rank 4
20210298 – Guidewire Clinical Survey Report (2020)	Angio-Seal™ VIP VCD	Rank 4
Test sugli animali post-commercializzazione	Iteration of Angio-Seal	MDCG 2020-06 Rank
Angio-Seal® Vascular Closure Device Biological Evaluation & Testing Report 17-3-121 (2018)	Angio-Seal™ VIP VCD Angio-Seal™ Evolution	Rank 12
Animal Study Report, Evolution CER044, Chronic Animal Study 93d, APS #90014675 (2008)	Angio-Seal™ VIP VCD Angio-Seal™ Evolution	Rank 11
Animal Study Report, Evolution CER056, Acute Animal Study at APS #90017642 (2008)	Angio-Seal™ VIP VCD Angio-Seal™ Evolution	Rank 11
Animal Study Report, Evolution CER056, Acute Animal Study at APS #90018838 (2008)	Angio-Seal™ VIP VCD Angio-Seal™ Evolution	Rank 11
Studi clinici post-commercializzazione pianificati/in corso	Iteration of Angio-Seal	MDCG 2020-06 Rank
The Angio-Seal Close: Post-market Clinical Investigation of the Angio-Seal™ VIP Vascular Closure Device (VCD): A Prospective, Multi-Center, Observational Study in Europe (T139E4)	Angio-Seal™ VIP VCD	Rank 3 (ongoing)

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1. Arunothayaraj S, Palmer S, Whitbourn R, et al. The Art of the Seal. J Am Coll Cardiol Interv. 2021 Jun, 14 (11) e117–e118. <https://doi.org/10.1016/j.jcin.2020.12.043>
2. E Kim, BG Sebastiao, A Lee, S Ande, J Shankar. Safety and effectiveness of vascular closure devices in interventional radiological procedures. Interventional neuroradiology: journal of peritherapeutic neuroradiology, surgical procedures and related neurosciences. 2022. :#pages#
3. F Frenzel, P Fries, R Shayesteh-Kheslat, A Buecker, A Massmann. Single Angio-Seal vascular closure device for transfemoral access exceeding 8F. Journal of cardiology. 2020. 76:#pages#
4. G Zenunaj, L Traina, P Acciarri, M Mucignat, S Scian, F Alesiani, R Serra, V Gasbarro. Superficial Femoral Artery Access for Infringuinal Antegrade Endovascular Interventions in the Hostile Groin: A Prospective Randomized Study. Annals of vascular surgery. 2022. :#pages#
5. Kreutz, Rolf P., et al. "Femoral artery closure devices vs manual compression during cardiac catheterization and percutaneous coronary intervention." Journal of the Society for Cardiovascular Angiography & Interventions (2022): 100370.
6. LE Jones, KH Yang, RW Feldtman, PV Uceda, CA Ferrara, JM Caruso, JL Richmond, SS Ahn. Safety and Efficacy of Arterial Closure Devices in an Office-Based Angiosuite. Annals of vascular surgery. 2018. 51:#pages#
7. N Ito, K Hishikari, H Hikita, A Takahashi, T Yonetsu, T Sasano. Images of angioscopy and intravascular ultrasound for dislodged Angio-Seal VIP vascular closure device. Cardiovascular intervention and therapeutics. 2022. 37:#pages#
8. M. Shukor, et al. "Managing Unintentional Iatrogenic Triple Lumen Central Venous Catheter (TLCVC) Arterial Puncture." The Interventionalist Journal 1.2 (2021): 12-18.
9. PS Dhillon, N McConachie, R Lenthall, S Nair, W Izzath. Use of 8Fr angio-seal for closure of femoral arteriotomy following use of 8Fr and 9Fr sheaths in patients undergoing mechanical thrombectomy for acute ischaemic stroke. Journal of clinical neuroscience : official journal of the Neurosurgical Society of Australasia. 2021. 88:#pages#

10. Shah, Aakash, et al. "Complete percutaneous decannulation from femoral venoarterial extracorporeal membrane oxygenation. "JTCVS techniques 6 (2021): 75-81.
11. Woźniak, Witold, and Piotr Ciostek. "Infectious arteritis of the femoral artery as a cause of recurrent hemorrhage after endovascular treatment of thrombosed popliteal artery aneurysm." *Acta Angiologica* 26.3 (2020): 113-117.
12. Su SF, Chang MY, Wu MS, Liao YC. Safety and efficacy of using vascular closure devices for hemostasis on sheath removal after a transfemoral artery percutaneous coronary intervention. *Jpn J Nurs Sci*. 2019 Apr;16(2):172-183. doi: 10.1111/jjns.12221. Epub 2018 Jul 25.

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Potvrđne studije pre plasiranja na tržište – originalna iteracija (HPCD)	Iteracija Angio-Seal	MDCG 2020-06 Rank
US Phase II European Phase II Early Ambulation in Diagnostic Angiography Patients	Angio-Seal Hemostatic Closure Devices (originalna iteracija)	Rank 4
Studije za odobrenje pre plasiranja na tržište (koje vode do Angio-Seal VIP VCD)	Iteracija Angio-Seal	MDCG 2020-06 Rank
Clinical Evaluation of the Angio-Seal Self-Tightening Slipknot (STS) Hemostatic Puncture Closure Device (IDE G990306)	Angio-Seal STS	Rank 2
Angio-Seal STS Early Ambulation and Discharge Study	Angio-Seal STS	Rank 2
Završene kliničke studije nakon plasiranja na tržište	Iteracija Angio-Seal	MDCG 2020-06 Rank
Clinical Study of Angio-Seal™ VIP Vascular Closure Device following Diagnostic and/or Interventional Endovascular Procedures (CP0901) (2009)	Angio-Seal™ VIP VCD	Rank 2
SEAL PM Post-Market Study of the St. Jude Medical Angio-Seal™ VIP Vascular Closure Device CV-13-008-EU-VC (2013)	Angio-Seal™ VIP VCD	Rank 2
Završene kliničke aktivnosti praćenja nakon plasiranja na tržište	Iteracija Angio-Seal	MDCG 2020-06 Rank
20210293 – Angio-Seal VIP Vascular Closure Device Clinical Survey Report (2020)	Angio-Seal™ VIP VCD	Rank 4
20210298 – Guidewire Clinical Survey Report (2020)	Angio-Seal™ VIP VCD	Rank 4
Testiranje na životinjama nakon plasiranja na tržište	Iteracija Angio-Seal	MDCG 2020-06 Rank
Angio-Seal® Vascular Closure Device Biological Evaluation & Testing Report 17-3-121 (2018)	Angio-Seal™ VIP VCD Angio-Seal™ Evolution	Rank 12
Animal Study Report, Evolution CER044, Chronic Animal Study 93d, APS #90014675 (2008)	Angio-Seal™ VIP VCD Angio-Seal™ Evolution	Rank 11
Animal Study Report, Evolution CER056, Acute Animal Study at APS #90017642 (2008)	Angio-Seal™ VIP VCD Angio-Seal™ Evolution	Rank 11
Animal Study Report, Evolution CER056, Acute Animal Study at APS #90018838 (2008)	Angio-Seal™ VIP VCD Angio-Seal™ Evolution	Rank 11
Planirane/tekuće kliničke studije nakon plasiranja na tržište	Iteracija Angio-Seal	MDCG 2020-06 Rank
The Angio-Seal Close: Post-market Clinical Investigation of the Angio-Seal™ VIP Vascular Closure Device (VCD): A Prospective, Multi-Center, Observational Study in Europe (T139E4)	Angio-Seal™ VIP VCD	Rank 3 (ongoing)

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1. Arunothayaraj S, Palmer S, Whitbourn R, et al. The Art of the Seal. J Am Coll Cardiol Interv. 2021 Jun, 14 (11) e117–e118. <https://doi.org/10.1016/j.jcin.2020.12.043>
2. E Kim, BG Sebastiao, A Lee, S Ande, J Shankar. Safety and effectiveness of vascular closure devices in interventional radiological procedures. Interventional neuroradiology: journal of peritherapeutic neuroradiology, surgical procedures and related neurosciences. 2022. :#pages#
3. F Frenzel, P Fries, R Shayesteh-Kheslat, A Buecker, A Massmann. Single Angio-Seal vascular closure device for transfemoral access exceeding 8F. Journal of cardiology. 2020. 76:#pages#
4. G Zenunaj, L Traina, P Acciarri, M Mucignat, S Scian, F Alesiani, R Serra, V Gasbarro. Superficial Femoral Artery Access for Infrainguinal Antegrade Endovascular Interventions in the Hostile Groin: A Prospective Randomized Study. Annals of vascular surgery. 2022. :#pages#
5. Kreutz, Rolf P., et al. "Femoral artery closure devices vs manual compression during cardiac catheterization and percutaneous coronary intervention." Journal of the Society for Cardiovascular Angiography & Interventions (2022): 100370.
6. LE Jones, KH Yang, RW Feldtman, PV Uceda, CA Ferrara, JM Caruso, JL Richmond, SS Ahn. Safety and Efficacy of Arterial Closure Devices in an Office-Based Angiosuite. Annals of vascular surgery. 2018. 51:#pages#
7. N Ito, K Hishikari, H Hikita, A Takahashi, T Yonetsu, T Sasano. Images of angiography and intravascular ultrasound for dislodged Angio-Seal VIP vascular closure device. Cardiovascular intervention and therapeutics. 2022. 37:#pages#
8. M. Shukor, et al. "Managing Unintentional Iatrogenic Triple Lumen Central Venous Catheter (TLCVC) Arterial Puncture."The Interventionalist Journal 1.2 (2021): 12-18.
9. PS Dhillon, N McConachie, R Lenthall, S Nair, W Izzath. Use of 8Fr angio-seal for closure of femoral arteriotomy following use of 8Fr and 9Fr sheaths in patients undergoing mechanical thrombectomy for acute ischaemic stroke. Journal of clinical neuroscience : official journal of the Neurosurgical Society of Australasia. 2021. 88:#pages#
10. Shah, Aakash, et al. "Complete percutaneous decannulation from femoral venoarterial extracorporeal membrane oxygenation. "JTCVS techniques 6 (2021): 75-81.

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1. Arunothayaraj S, Palmer S, Whitbourn R, et al. The Art of the Seal. J Am Coll Cardiol Interv. 2021 Jun, 14 (11) e117–e118. <https://doi.org/10.1016/j.jcin.2020.12.043>
2. E Kim, BG Sebastiao, A Lee, S Ande, J Shankar. Safety and effectiveness of vascular closure devices in interventional radiological procedures. Interventional neuroradiology: journal of peritherapeutic neuroradiology, surgical procedures and related neurosciences. 2022. :#pages#
3. F Frenzel, P Fries, R Shayesteh-Kheslat, A Buecker, A Massmann. Single Angio-Seal vascular closure device for transfemoral access exceeding 8F. Journal of cardiology. 2020. 76:#pages#
4. G Zenunaj, L Traina, P Acciarri, M Mucignat, S Scian, F Alesiani, R Serra, V Gasbarro. Superficial Femoral Artery Access for Infrainguinal Antegrade Endovascular Interventions in the Hostile Groin: A Prospective Randomized Study. Annals of vascular surgery. 2022. :#pages#
5. Kreutz, Rolf P., et al. "Femoral artery closure devices vs manual compression during cardiac catheterization and percutaneous coronary intervention." Journal of the Society for Cardiovascular Angiography & Interventions (2022): 100370.
6. LE Jones, KH Yang, RW Feldtman, PV Uceda, CA Ferrara, JM Caruso, JL Richmond, SS Ahn. Safety and Efficacy of Arterial Closure Devices in an Office-Based Angiosuite. Annals of vascular surgery. 2018. 51:#pages#
7. N Ito, K Hishikari, H Hikita, A Takahashi, T Yonetsu, T Sasano. Images of angioscopy and intravascular ultrasound for dislodged Angio-Seal VIP vascular closure device. Cardiovascular intervention and therapeutics. 2022. 37:#pages#
8. M. Shukor, et al. "Managing Unintentional Iatrogenic Triple Lumen Central Venous Catheter (TLCVC) Arterial Puncture."The Interventionalist Journal 1.2 (2021): 12-18.
9. PS Dhillon, N McConachie, R Lenthall, S Nair, W Izzath. Use of 8Fr angio-seal for closure of femoral arteriotomy following use of 8Fr and 9Fr sheaths in patients undergoing mechanical thrombectomy for acute ischaemic stroke. Journal of clinical neuroscience : official journal of the Neurosurgical Society of Australasia. 2021. 88:#pages#
10. Shah, Aakash, et al. "Complete percutaneous decannulation from femoral venoarterial extracorporeal membrane oxygenation."JTCVS techniques 6 (2021): 75-81.

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2.0 BIBLIOGRAFIE

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1. Arunothayaraj S, Palmer S, Whitbourn R, et al. The Art of the Seal. J Am Coll Cardiol Interv. 2021 Jun, 14 (11) e117–e118. <https://doi.org/10.1016/j.jcin.2020.12.043>
2. E Kim, BG Sebastiao, A Lee, S Ande, J Shankar. Safety and effectiveness of vascular closure devices in interventional radiological procedures. Interventional neuroradiology: journal of peritherapeutic neuroradiology, surgical procedures and related neurosciences. 2022. :#pages#
3. F Frenzel, P Fries, R Shayesteh-Kheslat, A Buecker, A Massmann. Single Angio-Seal vascular closure device for transfemoral access exceeding 8F. Journal of cardiology. 2020. 76:#pages#
4. G Zenunaj, L Traina, P Acciarri, M Mucignat, S Scian, F Alesiani, R Serra, V Gasbarro. Superficial Femoral Artery Access for Infrainguinal Antegrade Endovascular Interventions in the Hostile Groin: A Prospective Randomized Study. Annals of vascular surgery. 2022. :#pages#
5. Kreutz, Rolf P., et al. "Femoral artery closure devices vs manual compression during cardiac catheterization and percutaneous coronary intervention." Journal of the Society for Cardiovascular Angiography & Interventions (2022): 100370.
6. LE Jones, KH Yang, RW Feldtman, PV Uceda, CA Ferrara, JM Caruso, JL Richmond, SS Ahn. Safety and Efficacy of Arterial Closure Devices in an Office-Based Angiosuite. Annals of vascular surgery. 2018. 51:#pages#
7. N Ito, K Hishikari, H Hikita, A Takahashi, T Yonetsu, T Sasano. Images of angioscopy and intravascular ultrasound for dislodged Angio-Seal VIP vascular closure device. Cardiovascular intervention and therapeutics. 2022. 37:#pages#
8. M. Shukor, et al. "Managing Unintentional Iatrogenic Triple Lumen Central Venous Catheter (TLCVC) Arterial Puncture." The Interventionalist Journal 1.2 (2021): 12-18.
9. PS Dhillon, N McConachie, R Lenthall, S Nair, W Izzath. Use of 8Fr angio-seal for closure of femoral arteriotomy following use of 8Fr and 9Fr sheaths in patients undergoing mechanical thrombectomy for acute ischaemic stroke. Journal of clinical neuroscience : official journal of the Neurosurgical Society of Australasia. 2021. 88:#pages#
10. Shah, Aakash, et al. "Complete percutaneous decannulation from femoral venoarterial extracorporeal membrane oxygenation. "JTCVS techniques 6 (2021): 75-81.

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1. Arunothayaraj S, Palmer S, Whitbourn R, et al. The Art of the Seal. J Am Coll Cardiol Interv. 2021 Jun, 14 (11) e117–e118. <https://doi.org/10.1016/j.jcin.2020.12.043>
2. E Kim, BG Sebastiao, A Lee, S Ande, J Shankar. Safety and effectiveness of vascular closure devices in interventional radiological procedures. Interventional neuroradiology: journal of peritherapeutic neuroradiology, surgical procedures and related neurosciences. 2022. :#pages#
3. F Frenzel, P Fries, R Shayesteh-Kheslat, A Buecker, A Massmann. Single Angio-Seal vascular closure device for transfemoral access exceeding 8F. Journal of cardiology. 2020. 76:#pages#
4. G Zenunaj, L Traina, P Acciarri, M Mucignat, S Scian, F Alesiani, R Serra, V Gasbarro. Superficial Femoral Artery Access for Infrainguinal Antegrade Endovascular Interventions in the Hostile Groin: A Prospective Randomized Study. Annals of vascular surgery. 2022. :#pages#
5. Kreutz, Rolf P., et al. "Femoral artery closure devices vs manual compression during cardiac catheterization and percutaneous coronary intervention." Journal of the Society for Cardiovascular Angiography & Interventions (2022): 100370.
6. LE Jones, KH Yang, RW Feldtman, PV Uceda, CA Ferrara, JM Caruso, JL Richmond, SS Ahn. Safety and Efficacy of Arterial Closure Devices in an Office-Based Angiosuite. Annals of vascular surgery. 2018. 51:#pages#
7. N Ito, K Hishikari, H Hikita, A Takahashi, T Yonetsu, T Sasano. Images of angioscopy and intravascular ultrasound for dislodged Angio-Seal VIP vascular closure device. Cardiovascular intervention and therapeutics. 2022. 37:#pages#
8. M. Shukor, et al. "Managing Unintentional Iatrogenic Triple Lumen Central Venous Catheter (TLCVC) Arterial Puncture."The Interventionalist Journal 1.2 (2021): 12-18.
9. PS Dhillon, N McConachie, R Lenthall, S Nair, W Izzath. Use of 8Fr angio-seal for closure of femoral arteriotomy following use of 8Fr and 9Fr sheaths in patients undergoing mechanical thrombectomy for acute ischaemic stroke. Journal of clinical neuroscience : official journal of the Neurosurgical Society of Australasia. 2021. 88:#pages#
10. Shah, Aakash, et al. "Complete percutaneous decannulation from femoral venoarterial extracorporeal membrane oxygenation." JTCVS techniques 6 (2021): 75-81.

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11. Woźniak, Witold, and Piotr Ciostek. "Infectious arteritis of the femoral artery as a cause of recurrent hemorrhage after endovascular treatment of thrombosed popliteal artery aneurysm." *Acta Angiologica* 26.3 (2020): 113-117.
12. Su SF, Chang MY, Wu MS, Liao YC. Safety and efficacy of using vascular closure devices for hemostasis on sheath removal after a transfemoral artery percutaneous coronary intervention. *Jpn J Nurs Sci.* 2019 Apr;16(2):172-183. doi: 10.1111/jjns.12221. Epub 2018 Jul 25.

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İnsanlarda İlk Pazarlama Öncesi Angio-Seal - orijinal model (HPCD) Keşif Amaçlı Klinik Çalışma (8 French)	Iteration of Angio-Seal	MDCG 2020-06 Rank
Kensey Nash Hemostatic Puncture Closure Device (HPCD™) Phase I (G910028) (1992) 1. US Phase 1 2. US Phase 1 Extended 3. European Phase 1 4. European Phase 1 Extended	Angio-Seal Hemostatic Closure Devices (orijinal model)	Rank 2
İnsanlarda İlk Pazarlama Öncesi Angio-Seal - orijinal model (HPCD) Keşif Amaçlı Klinik Çalışma (6 French)	Iteration of Angio-Seal	MDCG 2020-06 Rank
6 French Pilot Study	Angio-Seal Hemostatic Closure Devices (orijinal model)	Rank 4
Pazarlama Öncesi Doğrulama Çalışmaları - orijinal model (HPCD)	Iteration of Angio-Seal	MDCG 2020-06 Rank
US Phase II European Phase II Early Ambulation in Diagnostic Angiography Patients	Angio-Seal Hemostatic Closure Devices (orijinal model)	Rank 4
Pazarlama Öncesi Onay Çalışmaları (Angio-Seal VIP VCD'ye Öncülük Eder)	Iteration of Angio-Seal	MDCG 2020-06 Rank
Clinical Evaluation of the Angio-Seal Self-Tightening Slipknot (STS) Hemostatic Puncture Closure Device (IDE G990306)	Angio-Seal STS	Rank 2
Angio-Seal STS Early Ambulation and Discharge Study	Angio-Seal STS	Rank 2

Tamamlanan Pazarlama Sonrası Klinik Çalışmalar	Iteration of Angio-Seal	MDCG 2020-06 Rank
Clinical Study of Angio-Seal™ VIP Vascular Closure Device following Diagnostic and/or Interventional Endovascular Procedures (CP0901) (2009)	Angio-Seal™ VIP VCD	Rank 2
SEAL PM Post-Market Study of the St. Jude Medical Angio-Seal™ VIP Vascular Closure Device CV-13-008-EU-VC (2013)	Angio-Seal™ VIP VCD	Rank 2
Tamamlanan Pazarlama Sonrası Klinik Takip Faaliyetleri	Iteration of Angio-Seal	MDCG 2020-06 Rank
20210293 – Angio-Seal VIP Vascular Closure Device Clinical Survey Report (2020)	Angio-Seal™ VIP VCD	Rank 4
20210298 – Guidewire Clinical Survey Report (2020)	Angio-Seal™ VIP VCD	Rank 4
Pazarlama Sonrası Hayvan Testi	Iteration of Angio-Seal	MDCG 2020-06 Rank
Angio-Seal® Vascular Closure Device Biological Evaluation & Testing Report 17-3-121 (2018)	Angio-Seal™ VIP VCD Angio-Seal™ Evolution	Rank 12
Animal Study Report, Evolution CER044, Chronic Animal Study 93d, APS #90014675 (2008)	Angio-Seal™ VIP VCD Angio-Seal™ Evolution	Rank 11
Animal Study Report, Evolution CER056, Acute Animal Study at APS #90017642 (2008)	Angio-Seal™ VIP VCD Angio-Seal™ Evolution	Rank 11
Animal Study Report, Evolution CER056, Acute Animal Study at APS #90018838 (2008)	Angio-Seal™ VIP VCD Angio-Seal™ Evolution	Rank 11
Planlanan / Devam Eden Pazarlama Sonrası Klinik Çalışmalar	Iteration of Angio-Seal	MDCG 2020-06 Rank
The Angio-Seal Close: Post-market Clinical Investigation of the Angio-Seal™ VIP Vascular Closure Device (VCD): A Prospective, Multi-Center, Observational Study in Europe (T139E4)	Angio-Seal™ VIP VCD	Rank 3 (ongoing)

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