ENGLISH – INSTRUCTIONS FOR USE.....	1
Angio-Seal™ VIP Vascular Closure Device	
DEUTSCH – GEBRAUCHSANLEITUNG	6
Angio-Seal™ VIP Vaskuläres Verschlusssystem	
FRANÇAIS – MODE D'EMPLOI.....	11
Système de fermeture vasculaire Angio-Seal™ VIP	
ITALIANO – ISTRUZIONI PER L'USO	16
Dispositivo di chiusura vascolare Angio-Seal™ VIP	
ESPAÑOL – INSTRUCCIONES DE USO	21
Dispositivo de cierre vascular Angio-Seal™ VIP	
ΕΛΛΗΝΙΚΑ – ΟΔΗΓΙΕΣ ΧΡΗΣΗΣ	26
Συσκευή αγγειακής αιμόστασης (σύγκλεισης) Angio-Seal™-VIP	
НА РУССКОМ ЯЗЫКЕ — ИНСТРУКЦИЯ ПО ЭКСПЛУАТАЦИИ	31
Устройство для закрытия места пункции сосуда Angio-Seal™ VIP	
MAGYAR – HASZNÁLATI UTASÍTÁS	36
Angio-Seal™ – VIP érelzáró eszköz	
PORTUGUÊS - INSTRUÇÕES DE UTILIZAÇÃO	41
Dispositivo de Encerramento Vascular Angio-Seal™ – VIP	
NEDERLANDS – GEBRUIKSAANWIJZING.....	46
Angio-Seal™ VIP Vaatafsluitsysteem	
SUOMI – KÄYTTÖOHJEET	51
Angio-Seal™ – VIP verisuonen sulkulaite	
SVENSKA – BRUKSANVISNING	56
Angio-Seal™ VIP Kärslutningsinstrument	
DANSK – BRUGERVEJLEDNING	61
Angio-Seal™ VIP – Udstyr til vaskulær forsegling	
NORSK – BRUKSANVISNING.....	66
Angio-Seal™ VIP vaskulær lukkeanordning	
POLSKI — INSTRUKCJA UŻYTKOWANIA.....	71
Urządzenie do zamykania naczyń Angio-Seal™ - VIP	
ČESKY – NÁVOD K POUŽITÍ	76
Vaskulární zastavovací zařízení VIP Angio-Seal™	
SLOVENSKY – NÁVOD NA POUŽITI	81
Cievne uzatváracie zariadenie Angio-Seal™ VIP	
EESTI KEEL – KASUTUSJUHEND	86
Angio-Seal™ VIP soonesulgur	
LIETUVIŲ K. – NAUDOJIMO INSTRUKCIJA	91
Kraujagyslių sandarinamasis prietaisas Angio-Seal™ - VIP	

TÜRKÇE – KULLANIM TALİMATLARI	96
Angio-Seal™ VIP Vasküler Kapatma Cihazı	
SRPSKI – UPUTSTVO ZA UPOTREBU	101
Angio-Seal™ - VIP uređaj za vaskularno zatvaranje	
한글 - 사용 설명서	106
Angio-Seal™ VIP 혈관 봉합 장치	
БЪЛГАРСКИ – ИНСТРУКЦИИ ЗА УПОТРЕБА	111
Изделие за съдово затваряне Angio-Seal™ VIP	
HRVATSKI – UPUTE ZA UPORABU	116
Uređaj za zatvaranje krvnih žila Angio-Seal™ VIP	
LATVIEŠU – LIETOŠANAS INSTRUKCIJA	121
Angio-Seal™ VIP asinsvadu slēgšanas ierīce	
LIMBA ROMÂNĂ - INSTRUCȚIUNI DE UTILIZARE	126
Dispozitiv de închidere vasculară Angio-Seal™ VIP	
SLOVENŠČINA – NAVODILA ZA UPORABO	131
Pripomoček za zapiranje žil VIP Angio-Seal™	
УКРАЇНСЬКА — ІНСТРУКЦІЯ З ВИКОРИСТАННЯ	136
Система для закриття місця пункції судин Angio-Seal™ VIP	

ENGLISH – INSTRUCTIONS FOR USE

Angio-Seal™ VIP Vascular Closure Device


TO ENSURE PROPER DEPLOYMENT AND USE OF THIS DEVICE AND TO PREVENT INJURY TO PATIENTS, READ ALL INFORMATION CONTAINED IN THESE INSTRUCTIONS FOR USE.

DEVICE DESCRIPTION

The Angio-Seal™ VIP Vascular Closure Device consists of the Angio-Seal device, an insertion sheath, an arteriotomy locator (modified dilator) and a guidewire with J-straightener. The Angio-Seal device is composed of a bovine-derived bioabsorbable collagen plug and a specially designed bioabsorbable polymer anchor that are connected by a bioabsorbable self-tightening suture (STS). The device seals and sandwiches the arteriotomy between its two primary members, the anchor and collagen plug. Hemostasis is achieved by the mechanical means of the anchor-arteriotomy-collagen sandwich. The device is contained in a delivery system that stores and then delivers the bioabsorbable components to the arterial puncture. The implanted components of the device are MRI Safe. The product is not made with natural rubber latex. The Angio-Seal VIP Vascular Closure Device functions as a system and does not require additional accessories beyond the guidewire with J-straightener included within the product's packaging for the device to function as intended.

Angio-Seal VIP Vascular Closure Device Components

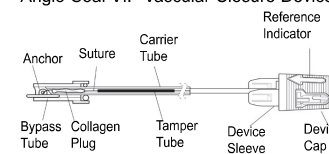PARTS OF THE ANGIO-SEAL VIP VASCULAR CLOSURE DEVICE

The Angio-Seal VIP Vascular Closure Device is supplied sterile in a poly bag. The bag includes an implant card and a tray containing the following components.

Unit	Part Name	Material
Bioabsorbable Unit *	Anchor	50:50 poly D,L lactide glycolide copolymer
	Suture	Polyglycolic acid Polycaprolate / polyglycolic acid copolymer coating
	Collagen Plug	Bovine hides
	Carrier Tube	Nylon or High-Density Polyethylene (HDPE)
Delivery System**	Tamper Tube	Polyurethane (6F) Polyethylene (8F)
	Insertion Sheath	Insertion Sheath Tubing Polyether Block Amine (PE-BAX) MDX/Hexane Solution
		Molding Acrylonitrile Butadiene Styrene (ABS)
		Sleeve Polyurethane (8F only)
		Hemostasis Seal Silicone/ Silicone Lubricant
		Insertion Sheath Cap Acrylonitrile Butadiene Styrene (ABS)
Arteriotomy Locator	Arteriotomy Locator	High Density Polyethylene (HDPE) Silicone lubricant
Accessories	Guidewire with a J-straightener 6Fr: 0.035" OD, 70 cm length guidewire 8Fr: 0.038" OD, 70 cm length guidewire	Guidewire: 304 Stainless Steel J-Straightener: Polypropylene

* Implantable portion of device

**Angio-Seal VIP Vascular Closure Device Delivery System consists of the delivery device containing the bioabsorbable unit and the insertion sheath.

The Angio-Seal VCD is sterilized via gamma sterilization and is a single-use medical device. The Angio-Seal device is designed in such a manner, that it cannot be re-used.

MODE OF ACTION

Hemostasis is achieved by sealing the arteriotomy with a bioabsorbable anchor and collagen plug. The anchor and collagen plug are held together (sandwiched) at the arteriotomy site by a bioabsorbable suture. The mechanism for deploying the collagen/anchor/suture uses a self-tightening slipknot to seal and sandwich the arteriotomy between the anchor and collagen plug. Hemostasis is achieved by the mechanical means of the anchor-arteriotomy-collagen sandwich. The ability of the implantable closure components to maintain compression on the arteriotomy has been tested at up to eight (8) hours.

ABSORBED MATERIALS

The Angio-Seal VIP Vascular Closure Device is composed of a bioabsorbable anchor, bioabsorbable hemostatic collagen plug and bioabsorbable suture. Hemostasis is achieved by properly forming a "sandwich" with the anchor located in the artery and the collagen plug outside connected through the suture. The anchor and suture will degrade by random hydrolytic cleavage of the ester function and the collagen plug will degrade through collagenase. The degradation products are metabolized and excreted in the urine or expired as carbon dioxide via the lungs. No accumulation effects have been observed in animal studies. The chemical composition, degradation products and absorption rate are described for each component in the following sections.

Suture

The chemical composition of the suture, braided polyglycolic acid with a coating of poly(glycolide-co-ε-caprolactone) 93:7 copolymer, the suture will degrade by random hydrolytic cleavage of the ester bonds forming the chemical entities listed in Table 1.

Table 1: Chemical composition of the polymer and hydrolysis products.

Component	Chemical Composition	Hydrolysis Product	Percent-age	Reaction
Suture	Polyglycolic acid braided suture	Glycolic Acid	100%	Monomer

Suture Coating	Glycolic acid	Glycolic Acid	93%	Monomer
	ε-Caprolactone	6-hydroxycaproic acid	7%	Monomer

Anchor

The anchor is molded from poly-(DL-lactide-co-glycolide) and will degrade by hydrolytic cleavage of ester bonds to form glycolic acid and lactic acid (Table 2).

Table 2: Chemical composition of the copolymer and hydrolysis product

Chemical Composition	Hydrolysis Product	Percentage	Reaction
Lactide, copolymer	Lactic acid	45-55%	Monomer
Glycolide, copolymer	Glycolic acid	45-55%	Monomer
Lactide, residual monomer	Lactic acid	0-1.5%	N/A
Residual monomer, Glycolide	Glycolic acid	<0.2%	N/A
Residual Tin	N/A	<100 ppm	Catalyst

Collagen Type 1, Bovine

The collagen plug is type I collagen sourced from bovine hide. The collagen material will be metabolized into its constituent amino acids (Table 3).

Table 3: Composition and metabolic end point

Chemical Composition	Metabolic Degradation	Percentage	Reaction
Collagen I, Bovine Hide	Amino acids	100%	N/A

Rate of Absorption:

Within ninety days (90) the anchor and suture are fully degraded via hydrolysis and the collagen is integrated into the local cellular matrix and ultimately metabolized into its constituent amino acids. The degradation products are absorbed by the surrounding tissue and excreted from the body in the urine and expiration via lungs.

Amount Absorbed: The bioabsorbable components are completely absorbed into the body, one hundred percent (100%).

PERFORMANCE CHARACTERISTICS

- Angio-Seal VIP Vascular Closure Device allows for closure of a common femoral artery having lumen diameter > 4mm, after percutaneous arterial access with a procedural sheath.
- Angio-Seal VIP Vascular Closure Device creates a seal by placing a bioabsorbable collagen plug and bioabsorbable anchor above and below the arteriotomy respectively, that are held together by a bioabsorbable suture. All implantable components resorb within 90 days.
- Hemostasis is achieved from mechanical means.
- Bovine derived collagen plug consists of a 9-Hole Weave V-Twist folding pattern to cover the arteriotomy.

CLINICAL BENEFITS

Angio-Seal VIP Vascular Closure Device provides the following clinical benefits relative to vascular closure after common femoral arterial puncture:

- Reduction in time to hemostasis¹
- Reduction in time to ambulation¹

1. Martin JL, Pratsos A, Magargee E, Mayhew K, Pensyl C, Nunn M, Day F, Shapiro T. A randomized trial comparing compression, Perclose Proglide and Angio-Seal VIP for arterial closure following percutaneous coronary intervention: the CAP trial. Catheter Cardiovasc Interv. 2008 Jan 1;71(1):1-5.

INTENDED USE

The Anglo-Seal Vascular Closure Device is a bioabsorbable implantable medical device intended for use in closing and reducing the time to hemostasis of common femoral arterial punctures following arterial access procedures.

INTENDED USERS

Intended users are physicians with training qualifying them to perform arterial access and closure for endovascular procedures through the common femoral artery and have participated in a Terumo Medical Corporation Anglo-Seal physician instruction program.

INDICATIONS FOR USE

The Anglo-Seal Vascular Closure Device is indicated for use in closing and reducing time to hemostasis at the common femoral arterial puncture site in patients who have undergone diagnostic angiography procedures or interventional procedures using an 8 French or smaller procedural sheath for the 8Fr Anglo-Seal device and a 6 French or smaller procedural sheath for the 6Fr Anglo-Seal device.

The Anglo-Seal Vascular Closure Device is also indicated for use to allow patients who have undergone diagnostic angiography and/or interventional procedures to ambulate safely as soon as possible after sheath removal and device placement.

CONTRAINDICATIONS

There are no contraindications to the use of this device. Attention is drawn to the warnings and precautions.

TARGET PATIENT POPULATION

The safety and effectiveness of the Anglo-Seal VIP Vascular Closure Device has been established in patients 18 years of age and older who have undergone an endovascular procedure via common femoral arterial access.

SPECIAL PATIENT POPULATIONS

The safety and effectiveness of the Anglo-Seal VIP Vascular Closure Device has not been established in the following patient populations:

- Patients who have known allergies to beef products, collagen, polyglycolic or polylactic acid polymers, and/or caprolactone.
- Patients with pre-existing autoimmune disease.
- Patients undergoing therapeutic thrombolysis.
- Patients with a vascular graft or stent at puncture site.
- Patients with uncontrolled hypertension (>180 mm Hg systolic).
- Patients with a high risk of bleeding such as a bleeding disorder, including thrombocytopenia (<100,000 platelet count), thrombasthenia, von Willebrand's disease, anemia (Hgb<10 mg/dl, Hct<30) or on anti-coagulation therapy.
- Patient population under 18 years old or others with small common femoral artery size (<4 mm in diameter). Small common femoral artery size may prevent the Anglo-Seal anchor from deploying properly in these patients.
- Patients who are pregnant or lactating.
- Patients with Ankle-Brachial Index (ABI) of <0.9.
- Patients with uncontrolled diabetes.
- Patients with evidence of a preexisting hematoma, arteriovenous fistula, or pseudoaneurysm at the access site prior to start of common femoral artery closure procedure.
- Patients experiencing cardiogenic shock (hemodynamic instability requiring intravenous medications or mechanical support) during or immediately post-catheterization.
- Patients with renal and/or hepatic impairment.
- Patients with other relevant co-morbidity.
- Patients with specific racial and/or ethnic origins.
- Patients with myocardial infarction within 72 hours.

WARNINGS

- Do not use if the temperature indicator dot on package has changed from light gray to dark gray or black.
- Do not use if the package is damaged or any portion of the package has been previously opened.
- Do not use if the items in the package appear damaged or defective in any way.
- Do not use where bacterial contamination of the procedural sheath or surrounding tissues may have occurred as this may result in an infection.
- Do not use if the procedural sheath has been placed through the superficial femoral artery and into the profunda femoris as this may result in collagen deposition into the superficial femoral artery. This may reduce blood flow through the vessel leading to symptoms of distal arterial insufficiency.
- Do not use if the puncture site is at or distal to the bifurcation of the superficial femoral and profunda femoris artery, as this may result in 1) the anchor catching on the bifurcation or being positioned incorrectly, and/or 2) collagen deposition into the vessel. These events may reduce blood flow through the vessel leading to symptoms of distal arterial insufficiency.
- Do not use if the puncture site is proximal to the inguinal ligament as this may result in a retroperitoneal bleeding/hematoma.
- Do not use if stenosis and/or significant plaque is present in the vicinity of the common femoral arterial puncture site.
- Avoid use in patients with an active systemic infection or evidence of an active infection at the planned procedural site.
- Do not use if patient has known allergies to the materials of the Anglo-Seal VIP Vascular Closure Device (defined in Device Description section above) as this may result in an allergic reaction.
- Keep away from sunlight and artificial light until planned use.

PRECAUTIONS

- The Anglo-Seal VIP Vascular Closure Device should only be used in patients who have undergone diagnostic angiography procedures or interventional procedures using 6 Fr Anglo-Seal with a 6 Fr or smaller procedural sheath and the 8Fr Anglo-Seal with an 8 Fr or smaller procedural sheath.
- Use a single wall puncture technique. Do not puncture the posterior wall of the artery.
- If a patient has had a procedural sheath left in place for longer than 8 hours, consideration should be given to the use of prophylactic antibiotics before insertion of the Anglo-Seal device.
- The Anglo-Seal device should be used within one hour of opening the foil pouch. The biodegradable components will begin to deteriorate upon exposure to ambient conditions.
- Observe sterile technique at all times when using the Anglo-Seal VIP Vascular Closure Device.
- The Anglo-Seal VIP Vascular Closure Device is for single use only and should not be reused or resterilized in any manner. The Anglo-Seal VIP Vascular Closure device is designed in such a manner, that it cannot be re-used.
- The Anglo-Seal device must be inserted through the insertion sheath provided in the package. Do not substitute any other sheath.
- Use only the arteriotomy locator provided in the package to locate the puncture in the arterial wall.
- The guidewire provided in the packaging is only to be used for the Anglo-Seal VIP Vascular Closure Device procedure.
- Follow physician orders regarding patient ambulation and discharge.
- If the Anglo-Seal device does not anchor in the artery due to improper orientation of the anchor or patient vascular anatomy, the bioabsorbable components and delivery system should be withdrawn from the patient. Hemostasis can then be achieved by applying manual pressure.
- If repuncture at the same location of previous Anglo-Seal VIP Vascular Closure Device use is necessary in ≤ 90 days, re-entry 1 cm proximal to the previous access site can be performed. If any other vascular

closure devices were used in the target artery then refer to the respective patient implant card for any precautions. Before considering Anglo-Seal use, a common femoral angiogram of the site is indicated.

- If patients have clinically significant peripheral vascular disease, the Anglo-Seal VIP Vascular Closure Device can be deployed safely in patient arteries > 5 mm diameter when there is found to be no luminal narrowing of 40% or greater within 5 mm of the puncture site.¹
- Dispose of contaminated device, components, and packaging materials utilizing standard hospital procedures and universal precautions for biohazardous waste.

¹ Abando, A., Hood D., Weaver, F., and S. Katz. 2004. "The use of the Anglo-Seal device for femoral artery closure". *J. Vasc Surg* 40:287-90.

PRECAUTIONS AT TIME OF DISCHARGE

Before considering discharge per facility discharge protocol, assess the patient for the following clinical conditions:

- Bleeding and/or hematoma at the closure site.
- Pain while walking.
- Signs of infection at the closure site.

ADVERSE EVENTS, HARMS, AND COMPLICATIONS

The following adverse reactions or conditions may also be associated with one or more Anglo-Seal VIP Vascular Closure Device components (i.e., collagen, synthetic bioabsorbable suture, silicone-based coating, and/or synthetic bioabsorbable polymer):

- Allergic Reaction
- Aneurysm
- AV Fistula
- Blood Loss/Bleeding
- Bovine Spongiform Encephalopathy (BSE) / Transmissible Spongiform Encephalopathy (TSE) and other Zoonotic diseases
- Death
- Device Failure
- Ecchymosis
- Embolism
- Foreign Body Reaction
- Hematoma
- Hemorrhage
- Infection
- Inflammatory Reaction
- Numbness
- Pain
- Patient Discomfort
- Procedure delay
- Pseudoaneurysm
- Retroperitoneal Bleed
- Sepsis
- Thromboembolism
- Thrombosis
- Vessel occlusion/Lower limb ischemia
- Vessel Perforation
- Vessel Tissue Dissection/Laceration

ANGIO-SEAL VIP VASCULAR CLOSURE DEVICE INSERTION PROCEDURE

The medical techniques and procedures described in these Instructions for Use do not represent ALL medically acceptable protocols, nor are they intended as a substitute for the clinician's experience and judgment in treating any specific patient.

The Anglo-Seal procedure is composed of three stages:

- A. Locate the Artery
- B. Set the Anchor
- C. Seal the Puncture

A. Locate the Artery

1. Assess the puncture site location and evaluate


the common femoral artery characteristics prior to placing the Anglo-Seal™ device by injecting contrast medium through the procedural sheath followed by an angiogram.

2. Using sterile technique, remove the Anglo-Seal device contents from the foil package by opening the package at the arrow indicator, taking care to pull the foil seal apart completely before removing the Anglo-Seal device.

NOTE: The Anglo-Seal device must be used within one hour after opening the foil pouch due to the moisture-sensitive nature of the product.


3. Insert the arteriotomy locator into the Anglo-Seal insertion sheath (Figure 1), making sure the two pieces snap together securely. To ensure proper orientation of the arteriotomy locator with the insertion sheath, the hub of the locator and the insertion sheath cap fit together only in the correct position. The reference indicator on the locator hub must align with the reference indicator on the insertion sheath cap.

Figure 1

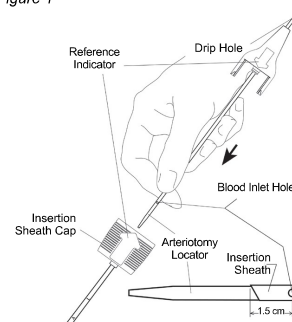4. Insert the Anglo-Seal guidewire into the procedural sheath that is currently in the patient. If the procedural sheath is smaller than the Anglo-Seal insertion sheath, it is advisable to ensure that the skin incision is of sufficient size to accommodate the Anglo-Seal insertion sheath.
5. Remove the procedural sheath, leaving the guidewire in place to maintain vascular access.
6. Thread the Anglo-Seal arteriotomy locator/insertion sheath assembly over the guidewire; the drip hole (located above the arteriotomy locator hub) will be oriented down and away so the flow of blood can be observed. Make sure the reference indicator on the insertion sheath is facing up, insert the assembly into the puncture tract. When the tip of the insertion sheath is about 1.5 cm into the artery blood will begin to flow from the drip hole in the locator (Figure 2).


Figure 2

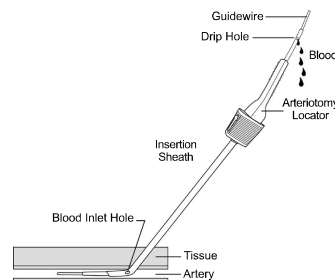NOTE: When resistance is encountered at the anterior vessel wall during over-the-guidewire advancement of the Anglo-Seal locator/insertion sheath assembly, rotate the assembly 90 degrees so the reference indicator is facing away from the user. This places the beveled tip of the insertion sheath perpendicular to the


- anterior vessel wall.
- Slowly withdraw the arteriotomy locator/insertion sheath assembly until blood slows or stops flowing from the drip hole. This indicates that the distal locator holes of the Angio-Seal insertion sheath have just exited the artery (Figure 3).

Figure 3

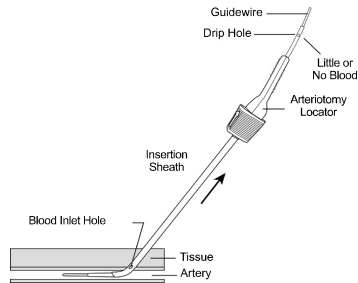
- From this point, advance the arteriotomy locator/insertion sheath assembly until blood begins to flow from the drip hole on the locator.
- NOTE: Over insertion of the arteriotomy locator/insertion sheath assembly into the artery, beyond 2 cm, may increase the chance of premature anchor hook up or interfere with the anchor's performance to achieve hemostasis.
- If blood flow does not resume, repeat Steps A-7 and A-8 until blood flows from the drip hole again upon advancement of the assembly into the artery.
- Holding the insertion sheath steady, without moving it into or out of the artery,
 - remove the arteriotomy locator and guidewire from the insertion sheath by flexing the arteriotomy locator upward at the insertion sheath hub (Figure 4).

Figure 4

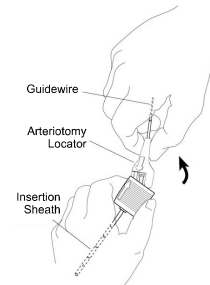- if necessary, rotate the insertion sheath so that the reference indicator on the insertion sheath cap is facing up (Figure 5).

Figure 5

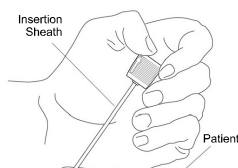WARNING: Under normal conditions, the Angio-Seal insertion sheath should not move into or out of the artery for the remainder of the Angio-Seal device deployment procedure. Using the insertion sheath markings as a guide, ensure that the insertion sheath position has not changed. If re-advancement is necessary, the guidewire and arteriotomy locator must be inserted prior to advancing the Angio-Seal insertion sheath.

B. Set the Anchor

- Confirm that the device sleeve has remained in the rear holding position (Figure 6). Carefully grasp the Angio-Seal device at the bypass tube. Cradle the Angio-Seal carrier tube in the palm of the hand and, with the reference indicator facing up, slowly insert the bypass tube and carrier tube into the insertion sheath hemostatic valve (Figure 7).

NOTE: If the device is not in the rear holding position (the tabs are past the detent), discard the device and replace with a new device.

Figure 6



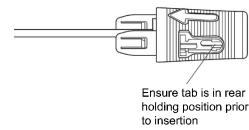
Figure 7

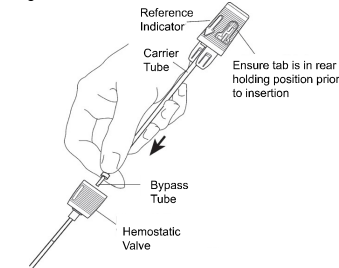- Confirm that the reference indicator on the insertion sheath is facing up. To ensure proper orientation of the Angio-Seal device with the insertion sheath, the insertion sheath cap and the device sleeve only fit together in the correct position. The reference indicator on the device cap should align with the reference indicator on the insertion sheath cap (Figure 8). Keeping the insertion sheath in place, carefully advance the Angio-Seal device in small increments until completely inserted into the insertion sheath. The insertion sheath cap and the device sleeve will snap together when properly fitted.

Figure 8

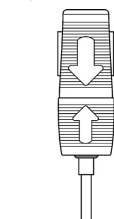
NOTE: If significant resistance to carrier tube advancement is encountered when insertion is almost complete, the anchor may be impinging on the posterior wall of the artery. DO NOT CONTINUE TO ATTEMPT TO ADVANCE. In this case, slight repositioning of the insertion sheath, either by reducing the angle of the insertion sheath with respect to the skin surface or by pulling the insertion sheath back by 1-2 mm, may permit normal deployment.

- With one hand continue to hold the insertion sheath cap steady to prevent movement of the insertion sheath into or out of the artery. With the other hand, grasp the device cap and slowly and carefully pull back. Slight resistance will be felt when the device cap is pulled out from rear holding position. Continue pulling on the device cap until resistance from the anchor catching on the distal tip of the insertion sheath is felt.

- To ensure the correct anchor position, confirm that the edge of the device cap falls within the colored bands on the device sleeve.

Figure 9 (Correct)


Resistance is felt as the carrier tube is withdrawn.

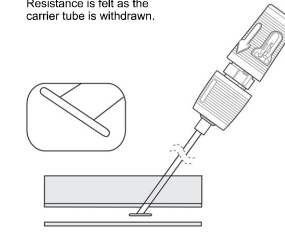- Maintain grip on the insertion sheath and pull the device handle straight back into the full rear locked position (Figure 10). Resistance will be felt as the device handle and sleeve snap lock. The device sleeve colored bands should now be completely visible.

NOTE: If the device sleeve separates from the insertion sheath while attempting full rear lock positioning, do not push the Angio-Seal™ device forward to reattach the insertion sheath cap. Complete sealing of the puncture following instructions in C-1.

Figure 10

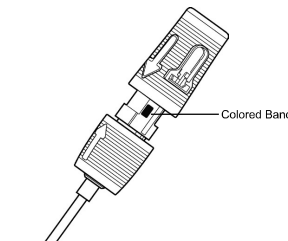
C. Seal the Puncture

- Once the anchor has been deployed correctly and the device cap has been locked into the rear position, slowly and carefully withdraw the device/insertion sheath assembly along the angle of the puncture tract to position the anchor against the vessel wall (Figure 11).


NOTE: Do not re-insert the device. Re-insertion of the device after partial deployment could cause collagen to be deposited in the artery.

Figure 11

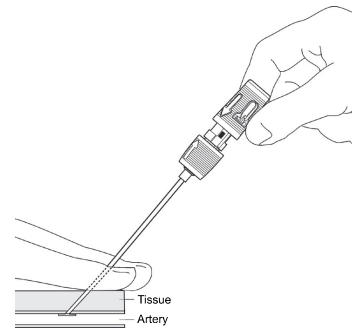- When the insertion sheath clears the skin, a tamper tube will appear (Figure 12).


Figure 12

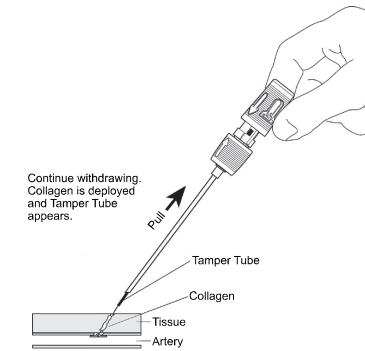WARNING: Failure to maintain tension on the suture while advancing the collagen could cause the collagen to enter the artery.

- Continue to withdraw the insertion sheath and device until the clear stop on the suture appears (Figure 13). Continue to pull until all of the suture has been deployed. The suture will then lock within the device cap where it is attached. Maintain tension on the suture.

Figure 13

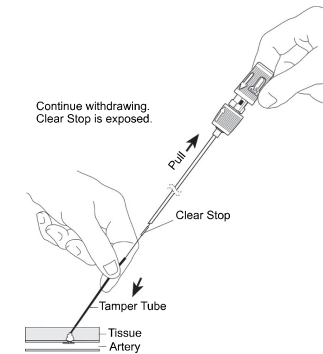- Continue to maintain tension on the suture during advancement of the knot and collagen with the tamper tube, following the puncture tract angle (Figure 14). A complete seal is indicated when resistance is felt, and hemostasis is achieved. As a guide, in most cases a compaction marker will be revealed. The essential indicators for a seal are resistance, hemostasis, and in most cases a compaction marker is revealed.

Figure 14


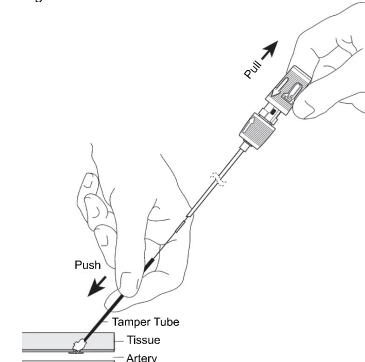
Figure 15

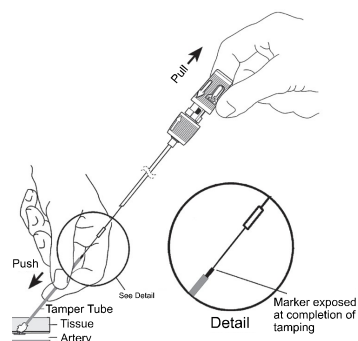
NOTE: The compaction marker is usually exposed. If not, repeat step C-4. Sufficient resistance will indicate completion of sealing and the compaction marker will usually be exposed (Figure 15).

NOTE: In some cases the compaction marker will not be fully exposed when resistance is felt. In such a case the clinician should check for hemostasis before completing the procedure.


NOTE: Once hemostasis is achieved, do not tamp to intentionally go beyond the distal end of the compaction marker (as shown in Figure 16) in order to prevent anchor deformation and/or collagen tearing.

Figure 16

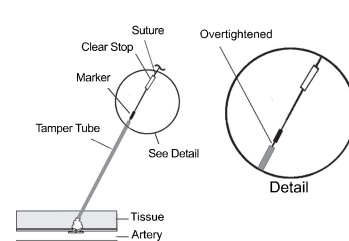
5. Cut the suture below the clear stop (Figure 17). Remove the tamper tube using a slight twirling upward motion.

Figure 17

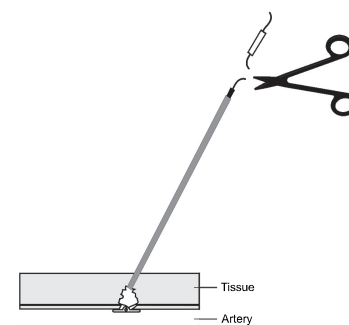6. Gently pull up on the suture. Push down on the skin using a sterile instrument. Cut the suture below the skin level, making sure to cut below the compaction marker (Figure 18).

Figure 18

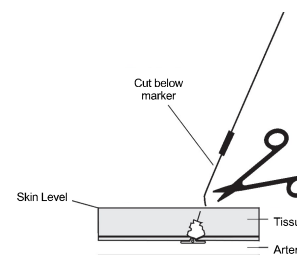The compaction marker is non-absorbable - it must be removed.

NOTE: Make sure the suture retracts below skin level to avoid infection.

NOTE: If seeping of blood occurs after placing the Angio-Seal™ device, or after removing the tamper tube, application of gentle digital pressure (one or two fingers) at the puncture site is usually sufficient to produce hemostasis. If manual pressure is necessary, monitor pedal pulses.

- Clean the puncture site with an antiseptic solution/ointment.
- Apply a sterile dressing to the puncture site so that it can be easily observed during recovery.
- While administering discharge instructions, provide the completed Patient Implant Card and Patient's Guide to the patient.

STORAGE AND PACKAGING

- Sterile and non-pyrogenic in unopened and undamaged package.
- This device is sterilized by Gamma Irradiation and is provided in a sterile state.
- The device may not be re-sterilized.
- Angio-Seal VIP Vascular Closure Device should be stored between 15°C to 25°C.
- Store the product in a manner that it is protected from sunlight and artificial light until planned use.

DEVICE DISPOSAL


The device should be disposed of using standard hospital procedures and universal precautions for biohazardous waste. All relevant laws and regulations governing such disposal should be followed.

SUMMARY OF SAFETY AND CLINICAL PERFORMANCE

The summary of safety and clinical performance (SSCP) is intended to enhance transparency by providing public access to an updated summary of clinical data and other information about the safety and clinical performance of the medical device. The SSCP is available in two versions, a version for healthcare professionals and a separate version for patients. The SSCP for the Angio-Seal VIP Vascular Closure Device medical device can be found at <https://www.terumomedical.com/ifu>

ANGIO-SEAL VIP VASCULAR CLOSURE DEVICE TROUBLESHOOTING

See table below for troubleshooting methods. Also refer to Adverse Events section for further information.

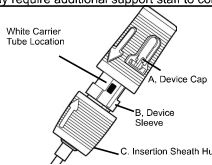Issues	Scenario	Possible Causes	Possible Actions	Preventive Actions
Collagen Exposure (Partial or Complete) / Thin Patients	Collagen is exposed above the skin after tamping.	Improper sealing technique. Narrow or fibrous tissue tract. Superficial common femoral artery.	Keep slight tension on the suture. Wet the collagen with saline. Attempt to advance the collagen under the skin with the tamper tube or hemostat (clamp). Tenting of the skin. Generate a lidocaine wheel at the puncture site. With collagen wet, attempt to advance collagen under the skin with tamper tube or hemostat. DO NOT apply vigorous tamping as this may result in anchor fracture. DO NOT cut off the excess collagen, as the suture woven through the collagen may be cut and the integrity of the anchor/collagen sandwich could be compromised. If a venous sheath is in place, remove the venous sheath to allow the anatomy to relax to its original position.	For a thin or patient with a narrow subcutaneous tract, make an incision on the skin surface at the puncture site before use. Dissect the subcutaneous tissue prior to insertion.
Obese Patients	The tamper tube may not be long enough to be exposed or grasped at the skin.	Deep common femoral artery.	Place fingers on either side of the suture, compress the surrounding tissue, and attempt to expose the tamper tube. If necessary, a sterile hemostat may be used to grasp the tamper tube so the collagen can be tamped adequately. May require additional support staff to compress the tissue adequately to expose tamper tube.	Attempt to perform a panniculus/pannus retraction to free the groin of being obscured.
Suture Lock-Up	Suture cannot be released from spool in the device cap. Requires force higher than normal between the collagen plug and carrier tube.	Mechanical problems with the spool in the main body.	 Cut through the carrier tube (slim white tube) and suture (between the base of the (A) device cap and the top of the (B) device sleeve and (C) insertion sheath hub assembly). Remove the device cap and insertion sheath assembly. Remove the carrier tube. Pull back on the suture. While maintaining tension on the suture, advance the tamper tube until the compaction marker is visible. Check for hemostasis.	N/A
Difficult Insertion of Insertion Sheath and Arteriotomy Locator Assembly	Insertion of the arteriotomy locator and insertion sheath is difficult.	Angio-Seal device is larger Fr size than procedural sheath. Fibrotic or narrow tissue tract. Thickened or sclerotic arterial wall. Improper insertion angle.	Enlarge tissue tract with hemostat (if indicated). Rotate the assembly 90° so the reference indicator is facing away from the user. This places the beveled tip of the insertion sheath perpendicular to the anterior vessel wall. Gently advance the insertion sheath tip through the tissue structures. Always rotate arteriotomy locator and insertion sheath assembly up again, so that insertion sheath markings and arrow are facing up, prior to the Angio-Seal device deployment. Use Seldinger technique angulation: ~45°. This may allow tip of insertion sheath to center in vessel lumen.	Ensure proper angle is used during access. Ensure bypass tube is completely seated within the insertion sheath valve. If patients have clinically significant peripheral vascular disease, the Angio-Seal device can be deployed in patient arteries >5mm diameter when there is no luminal narrowing of 40% or greater within 5mm of the puncture site.

Issues	Scenario	Possible Causes	Possible Actions	Preventive Actions
Compaction Marker Not Visible	Compaction marker is not visible after deployment.	Tamping angle is incorrect. Inadequate withdrawal of the insertion sheath and device. Patient is thin. Collagen is not moistened enough in tissue tract. See section above on Collagen Exposure. Patient is obese. See section above on obese patients.	Check angulation. Use Seldinger technique angulation: ~45°. Continue to maintain tension on the suture while continuing to advance the knot and collagen with the tamper tube.	Maintain angulation throughout procedure. Ensure clear stop is visible and all of the suture has been deployed. Maintain tension while advancing the tamper tube.
Difficult Insertion of carrier tube	Significant resistance to carrier tube advancement is encountered.	The insertion sheath may be bent. The anchor may be impinging on the posterior wall of the artery.	DO NOT CONTINUE TO ATTEMPT TO ADVANCE. Slight repositioning of the insertion sheath, either by reducing the angle of the insertion sheath with respect to the skin surface (less vertical) or by pulling the insertion sheath back by 1-2 mm. If resistance continues, remove device and apply manual or mechanical pressure per facility protocol.	Inspect the device for bends/kinks prior to insertion. Ensure tab is in rear holding position prior to insertion. Open the foil pouch/package completely while on a flat surface. Then remove the device to ensure bypass tube remains on device tip.
Collagen Deposition/Thrombosis	Collagen deposition (full or partial) of collagen plug into the artery. Thrombosis is suspected as patient develops cold leg and/or loss of distal pulse.	A) Stenosis or calcification is present near puncture site B) Puncture near the bifurcation or the anchor on the bifurcation C) Tension is not maintained on the suture while tamping. D) Premature anchor hook-up E) Deployment in a small caliber vessel F) Over-tamping past compaction marker	Confirm by duplex ultrasound. Request vascular surgery/surgeon consult for potential treatment to include, but not limited to: Thrombolysis, Percutaneous Thrombectomy or Surgical Intervention.	A) If patients have clinically significant peripheral vascular disease, the Angio-Seal device can be deployed in patient arteries >5mm diameter when there is no luminal narrowing of 40% or greater within 5mm of the puncture site B) Do not use if the puncture site is at or distal to the bifurcation C) Maintain tension on the suture during deployment D) If there are noticeable side branches within close proximity to the puncture site, take caution E) Ensure minimum vessel size is >4mm. Reference IFU for special patient populations F) Do not tamp past the compaction marker. Sufficient resistance will indicate completion of sealing. See Step C-4.
Insertion Sheath Movement	A) Insertion sheath has pulled out less than 1.5cm. B) Insertion sheath has pulled out greater than 1.5cm. C) Insertion sheath has been advanced into the vessel.	The insertion sheath is moved during removal of the arteriotomy locator and guidewire.	A) If re-advancement is necessary, the guidewire and arteriotomy locator must be re-inserted. Refer to Step A-6 through A-9. B) Remove device and apply manual or mechanical pressure to the puncture site per facility protocol. C) Pull back to previously noted position utilizing the insertion sheath markings.	Hold the insertion sheath steady without moving it into or out of the artery. Remove the arteriotomy locator and guidewire from the insertion sheath by flexing the arteriotomy locator upward at the insertion sheath hub. Do not pull the arteriotomy locator straight out of the insertion sheath.
Suture breakage	A) Above the tamper tube B) Inside the tamper tube C) Below the skin surface	The device was pulled back forcefully.	A) Grasp the remaining suture immediately and continue the procedure from Step C-4. B) Remove tamper tube and pull back on suture until resistance is felt. Observe for bleeding at the site. Continue the procedure from step C-5. If hemostasis has not been achieved, apply manual or mechanical pressure per facility protocol. C) If hemostasis has not been achieved, apply manual or mechanical pressure per facility protocol. Confirm if anchor was released by duplex ultrasound. If ischemic symptoms appear, treatment options include observation, percutaneous embolectomy of the anchor or fragments, or surgical intervention.	Avoid pulling back with excessive force. Avoid vigorous tamping.
Over-Insertion	Insertion sheath is advanced beyond locating arteriotomy which may result in premature anchor hook-up.	Insertion sheath is advanced too far into the vessel.	Reference IFU Steps A-6 through A-8.	Do not advance the arteriotomy locator and insertion sheath beyond 2cm.
Device Non-Deployment	Device pulls out with sheath upon withdrawal	Deployment in tissue tract due to improper arteriotomy location or fibrotic tissue tract Arteriotomy larger than anchor	Remove device and apply manual or mechanical pressure per facility protocol. Examine the device to ensure all bioabsorbable components have been withdrawn.	If while gaining access, difficult insertion or advancement of the procedural sheath is noted due to scar tissue or fibrotic tissue tract, deployment may be difficult.
Device Sleeve Separation	During the rear lock step, the device sleeve separates from its previously locked position within the insertion sheath cap.	Upon attempting to lock in rear lock position, operator rocks or wiggles device cap.	Grasp insertion sheath cap and withdraw device until tension is felt.	Do not rock or wiggle the device cap when pulling into the rear lock position.
No blood visible at the drip hole	No blood return during arteriotomy location step.	Arteriotomy locator and insertion sheath were not snapped together fully. Arteriotomy locator was inserted into insertion sheath backwards.	Maintain wire access and remove arteriotomy locator/insertion sheath assembly. Ensure components are snapped together properly and/or the reference indicators (arrows) are aligned properly per Step A-3. Reinsert over wire and continue from Step A-6.	When assembling arteriotomy locator and insertion sheath ensure arrows are aligned and that the arteriotomy locator snaps into the insertion sheath securely.
Component Breakage	The components may become brittle or break	Light Exposure Storage outside of labeled temperature range	Replace device	Review the temperature indicator on the sterile pouch and ensure that it is still light gray/silver. Do not store outside of required temperature requirements. Store devices in the original cardboard box free from light exposure. Do not store in a pyxis system. Do not store on a shelf exposed to artificial or sunlight
Bleeding/Blood Loss and Hematoma Post Deployment	Bleeding noticed at the closure site after intended positioning and deployment of the device.	Incomplete seal/insufficient tamping Required compressive sealing force not achieved.	Apply gentle digital or manual pressure per facility protocol to the puncture site. If manual pressure is necessary, monitor pedal pulses.	Tamp until the compaction marker is exposed. If compaction marker is not visible, sufficient resistance after tamping will ensure sealing.

IN CASE OF ADVERSE EVENTS, REPORT TO

Report adverse events to the legal manufacturer, EC Representative, and your EC Competent Authority of your member state in which the user and or patient is established.

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(C) 2021 Terumo Medical Corporation. All rights reserved. All brand names are trademarks or registered trademarks of their respective owners.
Terumo Australia Pty Ltd
Terumo Corporation
Terumo China (Hong Kong) Ltd:
MADE IN THE U.S.A.